

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification <sup>6</sup> : <b>C12Q 1/68</b>		A1	(11) International Publication Number: <b>WO 96/21741</b> (43) International Publication Date: 18 July 1996 (18.07.96)
(21) International Application Number: PCT/IB96/00026 (22) International Filing Date: 12 January 1996 (12.01.96) (30) Priority Data: 08/373,127 13 January 1995 (13.01.95) US 08/435,684 5 May 1995 (05.05.95) US (60) Parent Application or Grant (63) Related by Continuation US 08/435,684 (CON) Filed on 5 May 1995 (05.05.95) (71) Applicant (for all designated States except US): CIBA CORN- ING DIAGNOSTICS CORP. [US/US]; 63 North Street, Medfield, MA 02052 (US). (72) Inventors; and (75) Inventors/Applicants (for US only): SANDHU, Gurpreet, S. [IN/US]; 110 Seventh Avenue S.W., Rochester, MN 55902 (US). KLINE, Bruce, C. [US/US]; 2315 Kline Lane S.W., Rochester, MN 55902 (US). (74) Agent: MORGENSTERN, Arthur, S.; 63 North Street, Med- field, MA 02052 (US).			(81) Designated States: AL, AM, AT, AU, AZ, BB, BG, BR, BY, CA, CH, CN, CZ, DE, DK, EE, ES, FI, GB, GE, HU, IS, JP, KE, KG, KP, KR, KZ, LK, LR, LS, LT, LU, LV, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, TJ, TM, TR, TT, UA, UG, US, UZ, VN, ARIPO patent (KE, LS, MW, SD, SZ, UG), Eurasian patent (AZ, BY, KZ, RU, TJ, TM), European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, ML, MR, NE, SN, TD, TG).  Published With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.
(54) Title: NUCLEIC ACID PROBES FOR THE DETECTION AND IDENTIFICATION OF FUNGI			
<p>SEQ ID: 3 to SEQ ID: 23</p>			
(57) Abstract			
<p>Nucleic acid probes and primers are described for detecting fungi that cause disease in humans and animals, as well as spoilage of food and beverages. These probes can detect rRNA, rDNA or polymerase chain reaction products from a majority of fungi in clinical, environmental or food samples. Nucleic acid hybridization assay probes specific for <i>Acremonium</i> sp., <i>Aspergillus clavatus</i>, <i>Aspergillus flavus</i>, <i>Aspergillus fumigatus</i>, <i>Aspergillus glaucus</i>, <i>Aspergillus nidulans</i>, <i>Aspergillus niger</i>, <i>Aspergillus ochraceus</i>, <i>Aspergillus terreus</i>, <i>Aspergillus unguis</i>, <i>Aspergillus ustus</i>, <i>Beauveria</i> sp., <i>Bipolaris</i> sp., <i>Blastoschizomyces</i> sp., <i>Blastomyces dermatitidis</i>, <i>Candida albicans</i>, <i>Candida glabrata</i>, <i>Candida guilliermondii</i>, <i>Candida kefyr</i>, <i>Candida krusei</i>, <i>Candida lusitanae</i>, <i>Candida parapsilosis</i>, <i>Candida tropicalis</i>, <i>Chrysosporium</i> sp., <i>Cladosporium</i> sp., <i>Coccidioides immitis</i>, <i>Cryptococcus neoformans</i> var <i>gattii</i> serotype B, <i>Cryptococcus neoformans</i> serotype A, <i>Cryptococcus laurentii</i>, <i>Cryptococcus terreus</i>, <i>Curvularia</i> sp., <i>Fusarium</i> sp., <i>Filobasidium capsuligenum</i>, <i>Filobasidiella</i> (<i>Cryptococcus</i>) <i>neoformans</i> var <i>baicillispora</i> serotype C, <i>Filobasidiella</i> (<i>Cryptococcus</i>) <i>neoformans</i> var <i>neoformans</i> serotype D, <i>Filobasidium uniguttulatum</i>, <i>Geotrichum</i> sp., <i>Histoplasma capsulatum</i>, <i>Malbranchea</i> sp., <i>Mucor</i> sp., <i>Paecilomyces</i> sp., <i>Penicillium</i> species, <i>Pseudallescheria boydii</i>, <i>Rhizopus</i> sp., <i>Sporothrix schenckii</i>, <i>Scopulariopsis brevicaulis</i>, <i>Scopulariopsis brumptii</i>, <i>Saccharomyces cerevisiae</i>, and <i>Trichosporon beigellii</i> are also described.</p>			

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AM	Armenia	GB	United Kingdom	MW	Malawi
AT	Austria	GE	Georgia	MX	Mexico
AU	Australia	GN	Guinea	NE	Niger
BB	Barbados	GR	Greece	NL	Netherlands
BE	Belgium	HU	Hungary	NO	Norway
BF	Burkina Faso	IE	Ireland	NZ	New Zealand
BG	Bulgaria	IT	Italy	PL	Poland
BJ	Benin	JP	Japan	PT	Portugal
BR	Brazil	KE	Kenya	RO	Romania
BY	Belarus	KG	Kyrgyzstan	RU	Russian Federation
CA	Canada	KP	Democratic People's Republic of Korea	SD	Sudan
CF	Central African Republic	KR	Republic of Korea	SE	Sweden
CG	Congo	KZ	Kazakhstan	SG	Singapore
CH	Switzerland	LI	Liechtenstein	SI	Slovenia
CI	Côte d'Ivoire	LK	Sri Lanka	SK	Slovakia
CM	Cameroon	LR	Liberia	SN	Senegal
CN	China	LT	Lithuania	SZ	Swaziland
CS	Czechoslovakia	LU	Luxembourg	TD	Chad
CZ	Czech Republic	LV	Latvia	TG	Togo
DE	Germany	MC	Monaco	TJ	Tajikistan
DK	Denmark	MD	Republic of Moldova	TT	Trinidad and Tobago
EE	Estonia	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	UG	Uganda
FI	Finland	MN	Mongolia	US	United States of America
FR	France	MR	Mauritania	UZ	Uzbekistan
GA	Gabon			VN	Viet Nam

5

NUCLEIC ACID PROBES FOR THE DETECTION AND IDENTIFICATION  
OF FUNGI

10

FIELD OF INVENTION

The inventions described and claimed herein relate to the design and composition of two nucleic acid probes capable of detecting many different fungal organisms in clinical, food, environmental and other samples. The inventions described and claimed herein also relate to the design and composition of probes capable of specifically detecting and identifying *Acremonium* sp., *Aspergillus clavatus*, *Aspergillus flavus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*, *Aspergillus ustus*, *Beauveria* sp., *Bipolaris* sp., *Blastoschizomyces* sp., *Blastomyces dermatitidis*, *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* var *gattii* serotype B, *Cryptococcus neoformans* serotype A, *Cryptococcus laurentii*, *Cryptococcus terreus*, *Curvularia* sp., *Fusarium* sp., *Filobasidium capsuligenum*, *Filobasidiella* (*Cryptococcus*) *neoformans* var *bacillispora* serotype C, *Filobasidiella* (*Cryptococcus*) *neoformans* var *neoformans* serotype D, *Filobasidium uniguttulatum*, *Geotrichum* sp., *Histoplasma capsulatum*, *Malbranchea* sp., *Mucor* sp., *Paecilomyces* sp., *Penicillium* species, *Pseudallescheria boydii*, *Rhizopus* sp., *Sporothrix schenckii*, *Scopulariopsis*

*brevicaulis* sp., *Scopulariopsis brumpti*, *Saccharomyces cerevisiae*, and *Trichosporon beigelii* in clinical, food, environmental and other samples.

Fungi are eukaryotic microorganisms that are universally distributed. While  
5 in nature fungi play a major role in the decomposition of plant materials, they are also responsible for spoilage of food, beverage and pharmaceutical preparations. Out of an estimated 100,000 species of fungi described by mycologists, approximately 150 species are pathogenic to man and animals. The increasing incidence of AIDS and the development of newer treatments for hematologic  
10 malignancies and organ transplants has lead to an increase in the number of immunocompromised patients. These patients have a high risk of developing fungal infections, which if not rapidly diagnosed and treated are capable of causing death in a matter of days. The number of antifungal drugs is limited and their toxic side effects on the patient are much higher than that of comparable antibacterial therapy.  
15 A rapid diagnosis of fungal infection and start of treatment is critical in these patients. Books by Kwon-Chung and Bennett, along with Sarosi and Davies, provide an overview into the medical importance of fungi.

Fungal organisms are identified by morphology and nutritional  
20 characteristics. Fungi may take anywhere from two days to several weeks to grow in culture and often the same organism can take radically different forms depending on the growth conditions. This makes timely identification difficult even for the classically trained expert and impedes the treatment of patients where rapid identification of genus and species is of medical advantage.

25

The incidence and distribution of major pathogenic fungi varies by geographic location. *Aspergillus fumigatus*, *Blastomyces dermatitidis*, *Candida albicans*, *Coccidioides immitis*, *Cryptococcus neoformans*, *Histoplasma capsulatum*,

*Paracoccidioides brasiliensis*, *Pseudallescheria boydii* and *Sporothrix schenckii* represent some of the leading causes of mycotic infections.

*Aspergillus fumigatus* is among the top three causes of systemic fungal infection treated in hospitals. It usually affects patients with organ transplants, acute leukemias and burns and can be rapidly fatal if not diagnosed quickly. With over 150 species of *Aspergillus* present in the soil, air and water, accurate detection of *Aspergillus fumigatus* becomes extremely important. *Aspergillus clavatus*, *Aspergillus flavus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis* and *Aspergillus ustus* represent a majority of *Aspergillus* species seen in clinical specimens and their presence can cause diagnostic difficulties. *Aspergillus flavus*, *Aspergillus fumigatus* and *Aspergillus niger* have been linked with disease in humans, with *Aspergillus fumigatus* being the predominant pathogen in North America. A few immunologic tests exist for *Aspergillus fumigatus* but these have limited sensitivity and specificity. There are also reports of development of polymerase chain reaction based tests for *Aspergillus fumigatus* based on the amplification of the *Asp f1* antigen gene and a ribosomal intergenic spacer (Spreadbury et. al.). The Spreadbury technique is based on the PCR amplification of a 401 bp fragment spanning the large subunit rRNA/intergenic spacer region. This relies on a pair of primers to specifically amplify DNA from *Aspergillus fumigatus* only, and is of no utility in identifying other fungi.

*Blastomyces dermatitidis* is present in the soil, usually in bird droppings and animal feces. Infections often occur at construction sites and the ensuing lung infiltration and pneumonitis are usually fatal in immunocompromised patients. Diagnosis by culture may take weeks, and the organism is occasionally mistaken for other fungi. Existing immunological diagnostic tests are unreliable, and there is a need for rapid and reliable DNA based diagnostic tests. Similarly, *Histoplasma*

*capsulatum* exists in the soil and is known to have infected at least 20% of the population of North America. Most infections start in the lung and resolve spontaneously, but may occasionally spread to other organs. AIDS patients represent a growing number of cases of Histoplasmosis. Diagnosis is difficult as

5 immunological tests are often negative during the first 4-6 weeks of infection. *Coccidioides immitis* is found in abundance in the soil in Southwestern United States. Dust storms, farming, building construction, earthquakes and even hiking have been linked with outbreaks of disease. Lung infection followed by cavitation and disseminated miliary coccidioidomycosis are seen. Meningitis is usually lethal,

10 and as with other fungi, mortality is highest in debilitated hosts. Four serotypes of *Cryptococcus neoformans* cause disease in humans. These are *Cryptococcus neoformans* serotype A, *Cryptococcus neoformans* var *gatti* serotype B, *Filobasidiella (Cryptococcus) neoformans* var *bacillispora* serotype C and *Filobasidiella (Cryptococcus) neoformans* var. *neoformans* serotype D. The

15 incidence of this disease is growing rapidly, with up to 10% of HIV infected people developing cryptococcosis. DNA probes capable of detecting all 4 serotypes are required for the early diagnosis and treatment for life threatening infections like cryptococcal meningitis. A report by Stockman et. al. discusses commercial tests for *Histoplasma*, *Blastomyces*, *Coccidioides*, and *Cryptococcus* based on the 18S

20 rRNA (Gen-Probe, Inc., San Diego, CA). The authors report sensitivities ranging from 87.8 to 100% and a specificity of 100%. One drawback of these probes is that these are used on rRNA extracted from fungal cultures. As some fungi may require up to 3 weeks to grow in culture, this technique cannot be used to expedite diagnosis until a culture becomes available.

25

*Candida albicans* is one of the most common causes of fungal infection in humans. It is present in the respiratory, gastrointestinal and female genital tract of healthy individuals, and acts as an opportunistic pathogen in debilitated individuals on steroid or chemotherapy. Diabetes mellitus and indwelling catheters are other

predisposing causes. Immunocompromised hosts show rapid hematogenous spread of fungi. Morbidity and mortality in untreated cases is high. *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis* and *Candida tropicalis* are also known to cause disease in humans. DNA probes capable of identifying these individual species would eliminate the need for multiple blood cultures and lengthy biochemical speciation.

Recent advances in molecular techniques have led to the approach of microbe detection and identification based upon the DNA sequence of ribosomal genes. Commonly used detection techniques include either direct amplification of the ribosomal DNA (rDNA) genes by the polymerase chain reaction, or reverse transcription of the ribosomal RNA (rRNA) into complementary DNA (cDNA) followed by polymerase chain reaction amplification of the cDNA. Ribosomes are composites of unique rRNA and protein species that function in the translation of messenger RNA into protein. Evolutionary studies are consistent with the interpretation that all extant life has evolved from a single organism. Thus, all cellular organisms contain rRNA and these rRNAs are related by evolution. The evolutionary process is such that each species of organism appears to have unique regions of sequence in its ribosomal genes. The presence of these unique species specific regions allows one to design DNA probes that under conditions of hybridization will specifically bind to, and identify the polymerase chain reaction amplified DNA from only one species of fungus. For the purposes of this application, the word "primer" is used to mean a nucleotide sequence which can be extended by template-directed polymerization, and "probe" is used to mean a nucleotide sequence capable of detecting its complementary sequence by hybridization. Also, for the purpose of this application, the phrase "nucleotide sequence" is intended to include either DNA or RNA forms or modification thereof. Furthermore, those versed in the art will recognize that primer sequences can be used as probes and vice versa. The use of nucleic acid hybridization to detect

specific nucleic acid sequences of interest is also described by Kohne (U.S. Patent 4,851,330, 7/1989).

In prokaryotes and eukaryotes, ribosomal RNA and the corresponding  
5 rDNA genes are identified by the size of the RNA. The sizes are related in terms of sedimentation velocity or S values. Thus, for prokaryotes the values are 5S, 16S, and 23S; and for eukaryotes the values are 5S, 5.8S, 18S and 28S. Because all ribosomes perform the same function which is essential for cell viability, ribosomal sequences are largely conserved, yet certain regions of each ribosomal species are  
10 subject to more variation without consequence to function. It is these hypervariable regions that allow one to identify different species amongst members of the same genus. As noted in the references, there are several reports where 5S, 18S and the intergenic spacer between 5.8S and 28S rDNA have been used for the detection and identification of fungi (Holmes et. al., Hopfer et. al., Lott et. al., Maiwald et. al.,  
15 Makimura et. al., Mitchell et. al., Nakamura et. al.). Holmes et. al. describe a PCR test based on the co-amplification of the 5S rDNA and an adjacent nontranscribed spacer region. This identifies only *Candida albicans* and detects other *Candida* species without identifying individual organisms. Hopfer et. al. and Maiwald et. al. both use universal primers to amplify 18S rDNA from several fungi including  
20 *Candida* sp., *Aspergillus fumigatus*, *Cryptococcus neoformans* and *Trichosporon* sp. These amplicons are digested with restriction enzymes and the cut fragments are sized by gel electrophoresis. This restriction fragment length polymorphism pattern enables them to identify most but not all organisms. This technique can be used on amplified DNA from a pure fungal culture. As clinical samples such as sputum  
25 usually contain multiple fungal organisms, this technique has little utility in diagnosis as multiple overlapping fragments obtained from a mix of fungi would be nearly impossible to interpret. Lott et. al. use the 5.8S RNA and the internal transcribed spacer (ITS2) to identify and speciate *Candida albicans* and related *Candida* species. Makimura amplifies a 687 bp fragment from the 18S rDNA of 25



medically important fungi and uses these in the diagnosis of *Candida albicans* in clinical samples. Mitchell uses nested PCR to amplify 5.8S and internal transcribed spacer (ITS) to identify *Cryptococcus neoformans*. No subsequent testing is done to verify the identity of the amplified DNA. Nakamura et. al. use 18S primers to detect  
5 *Aspergillus fumigatus* infections of the lung. Most protocols given in these references can only be used to detect an extremely limited number of fungi from a clinical specimen. Hopfer et. al. and Maiwald et. al. can identify multiple organisms from pure cultures, but their utility for clinical specimens containing multiple fungal species is limited at best.

10

United States patents have been issued to Weisburg et. al. for probes developed for the detection of 18S small subunit ribosomal RNA sequences in fungi. These probes will detect fungi from many species, but cannot be used easily to identify any single species. United States patents have also been issued to  
15 Milliman for probes developed for the specific detection of the bacteria *Staphylococcus aureus* based on the 16S ribosomal sequences. Hogan et. al. (European Pat. App. 0,272,009) describe one fungal probe for 18S rRNA and three fungal probes for 28S rRNA sequences. Two of these 28S probes detect several different fungi while the third probe detects *Candida krusei* from a limited panel of  
20 10 fungi. None of the 28S probes described by Hogan et. al. is related to any of the probes described in our invention. All probes claimed in our invention can be mapped within the first 900 base pairs of a 28S gene. The probes described by Hogan et. al. are located further 3' on the 28S sequence, between base pairs 1000 and 2000 (these numbers are comparable to the primary sequence of *Saccharomyces*  
25 *cerevisiae* 28S rRNA gene. Genbank accession number: J01355). Leclerc et. al. have published reports analyzing the phylogenetic relationship between fungi based on partial DNA sequences of several fungal 28S genes sequenced by them. Some of the organisms claimed to have been sequenced by Leclerc are the same as some organisms sequenced by us. These are *Sporothrix schenckii*, *Pseudallescheria*

*boydii*, *Blastomyces dermatitidis*, *Histoplasma capsulatum* and *Chrysosporium* sp. Leclerc et. al. have not published any sequence data in their report, and to the best of our knowledge, they have not made these sequences publically available in the GenBank. The reverse-complement sequence of their sequencing primer 401  
5 (TCCCTTTCAA CAATTTACG) overlaps our SEQ ID NO: 1 (GTGAAATTGT TGAAAGGGAA) by 19 nucleotides and their sequencing primer 636 (GGTCCGTGTT TCAAGACGG) overlaps our SEQ ID NO: 2 (GACTCCTTGG TCCGTGTT) by 10 nucleotides. We are aware of no reports in the literature of  
10 development of species specific diagnostic probes.

As discussed above, most present techniques for the molecular detection of fungi rely on the use of highly specific primers for the PCR amplification of only one fungal species. Those that employ "Universal" primers for a PCR amplification  
15 of DNA from multiple organisms, use post-PCR amplicon identification techniques that are useful only on pure cultures of fungi. These are not be able to identify fungi from a clinical specimen containing multiple fungal organisms. Our first aim was to develop "Universal" primers for the 28S gene. These primers would be capable of amplifying in a PCR, 28S rDNA from most fungi. Our subsequent aim was to  
20 develop species specific probes for fungi of interest, that would be used to analyze our "Universal" 28S amplicon. These species specific probes would be able to detect the presence of fungi of interest even in situations containing mixed fungal species.

25 One aspect of this invention is to provide nucleic acid primers capable of detecting 28S sequences from DNA or RNA of most fungi. These would be used as "Universal" primers in a polymerase chain reaction to amplify 28S sequences from any fungus present in clinical, food, environmental or other samples. These "Universal" primers would also be used to sequence the amplified DNA. The

sequence obtained would be used to identify the fungus by comparing with a database of known fungal sequences.

A second aspect of this invention is to provide nucleic acid probes capable  
5 of detecting and identifying, by nucleic acid hybridization, the pathogens  
*Aspergillus fumigatus*, *Blastomyces dermatitidis*, *Candida albicans*, *Coccidioides immitis*, *Cryptococcus neoformans*, *Histoplasma capsulatum*, *Aspergillus flavus*,  
*Aspergillus glaucus*, *Aspergillus niger*, *Aspergillus terreus*, *Candida glabrata*,  
*Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitaniae*,  
10 *Candida parapsilosis*, *Candida tropicalis*, *Pseudallescheria boydii*, *Sporothrix schenckii* and other species by use of any of several different formats. Additionally,  
nucleotide sequence information is provided to identify these pathogens and other  
fungi by DNA sequence comparison (Figure 2) or by the construction of additional  
probes.

15

#### SUMMARY OF THE INVENTION

Nucleic acid probes and primers are described for detecting fungi that cause  
disease in humans and animals, as well as spoilage of food and beverages. These  
20 probes can detect rRNA, rDNA or polymerase chain reaction products from a  
majority of fungi in clinical, environmental or food samples. Nucleic acid  
hybridization assay probes specific for *Aspergillus fumigatus*, *Blastomyces dermatitidis*, *Candida albicans*, *Coccidioides immitis*, *Cryptococcus neoformans*,  
*Histoplasma capsulatum*, *Aspergillus flavus*, *Aspergillus glaucus*, *Aspergillus niger*,  
25 *Aspergillus terreus*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*,  
*Candida krusei*, *Candida lusitaniae*, *Candida parapsilosis*, *Candida tropicalis*,  
*Pseudallescheria boydii*, *Sporothrix schenckii* and other species (Table 1 and Figure  
2) are also described.

## BRIEF DESCRIPTION OF THE DRAWINGS

Figure 1 represents the relative position of the sequences described on the 28S subunit of fungi.

- 5        Figures 2A, B, C and D together represent the multiple sequence alignment for (SEQ ID NO: 24) through (SEQ ID NO: 74).

## DETAILS OF THE INVENTION

- 10        Our first objective was to develop nucleic acid primers for use in a polymerase chain reaction to amplify 28S genes from all fungi likely to be present in a clinical sample. This amplified DNA would then be amenable to probing with several different species specific probes. Each one of these species specific probes would, under conditions of hybridization, anneal to 28S ribosomal DNA from only
- 15        one species of fungus, thereby detecting and identifying the species of fungus present in the clinical sample. The 28S gene was selected as a target because it had regions that were conserved among fungi and these would provide potential annealing sites for "universal" fungal probes. The ribosomal 28S genes were also expected to have hypervariable regions that would be unique enough to provide
- 20        sites for species specific probes. The large rRNA gene is called the 23S rRNA gene in prokaryotes and 28S in eukaryotes. This designation is based on the length and therefore the sedimentation coefficient of these rRNA molecules. Fungal large subunit rRNAs vary in size among different organisms and are often referred to as being 25S, 26S or 28S. Since fungi are eukaryotes, and to maintain uniformity in
- 25        this application, we shall refer to fungal large subunit rRNA as 28S rRNA.

Published sequences from *Cryptococcus neoformans*, two *Candida albicans*, *Saccharomyces cerevisiae* and two *Schizosaccharomyces pombe* 28S genes are approximately 3.5 kilobases in length (Genbank accession numbers:

L14068, L28817, X70659, J01355, Z19136 & Z19578). These four sequences were aligned, and a region of sequence variability was found clustered between coordinates 200 and 700 from the 5' end of these genes. As an initial starting point, two nucleic acid primers P1 (ATCAATAAGC GGAGGAAAAG) and P2 (CTCTGGCTTC ACCCTATTC) (see figure 1), capable of hybridizing to all 4 of the above mentioned organisms and not to human 28S sequences (GenBank accession number: M11167), were designed and used under low stringency hybridization conditions in a polymerase chain reaction to amplify approximately 800 base pairs of DNA spanning this hypervariable region from the following 34 fungi that were obtained from the Mayo Clinic fungal collection: *Acremonium* sp., *Aspergillus clavatus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*, *Aspergillus ustus*, *Beauveria* sp., *Bipolaris* sp., *Blastomyces dermatitidis*, *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* serotype A, *Curvularia* sp., *Geotrichum* sp., *Histoplasma capsulatum*, *Mucor* sp., *Penicillium* sp., *Pseudallescheria boydii*, *Saccharomyces cerevisiae*, *Sporothrix schenkii* and *Trichosporon beigelii*.

20

DNA was extracted from the fungi listed above by the following method. A loopful of fungal culture was scraped off a culture plate using a sterile inoculation loop. The fungus was added one milliliter of sterile water in a 1.5 ml Sarsted (Newton, North Carolina) screw cap microcentrifuge tube. This tube was placed in a boiling water bath for 20 minutes in order to lyse the fungus and release DNA from the cells. Two microliters of this whole cell lysate was used in a PCR to amplify 28S rDNA. All PCR amplifications were carried out as hot-start reactions in a 50 ul reaction volume using Perkin-Elmer (Norwalk, CT) 0.5 ml thin-wall polypropylene tubes and a Perkin-Elmer thermal cycler. Reagents added to the tube

initially were 2.5 ul of 10X PCR buffer (100 mM tris pH 8.3, 500 mM KCl, 15 mM MgCl<sub>2</sub>), 5.0 ul of 50% glycerol/1 mM cresol red, 8.0 ul of dNTP mix (1.25 mM each of dATP, dGTP, dTTP and dCTP), 12 picomoles of each nucleic acid primer and sterile water to make up a volume of 25 ul. A wax bead (Ampliwax Gem-100, 5 Perkin-Elmer) was added and the tubes heated to 77°C for 1 minute and cooled to room temperature to form a wax barrier. 2.5 ul of 10X PCR buffer, 5.0 ul of 50% glycerol/1 mM cresol red, 0.2 ul Taq polymerase (AmpliTaq 5U/ul, Perkin-Elmer) and 15.3 ul of sterile water was added to the tube along with 2.0 ul of DNA from the fungal whole cell lysate described above. 50 cycles of thermal cycling was 10 carried out at 94°C - 30 sec, 40°C - 1 min, 72°C - 2 min. The amplified DNA was electrophoresed and purified from a low melt agarose gel by tris buffered phenol pH 8.0, phenol/chloroform/isoamyl alcohol (25:24:1 by vol.) and 3 ether extractions, followed by isopropanol precipitation and 70% ethanol wash.

15 We completely sequenced both strands of DNA amplified from the organisms listed above. All sequencing was carried out on an Applied Biosystems 373A sequencer. Every nucleotide in the sequences generated was verified and confirmed by examining the complementary nucleotide from the second strand sequence. We had now created a novel database consisting of nucleic acid 20 sequences spanning a variable region of the 28S rDNA from a diverse collection of medically important fungi.

While the complete sequences for *Candida albicans*, *Cryptococcus neoformans* and *Saccharomyces cerevisiae* 28S genes had previously been 25 published and deposited in GenBank, it was not obvious, nor had it been defined, whether any regions of sequence identity among these three organisms would also be conserved among all fungi of interest. DNA sequences from all the fungi in our novel 28S database had to be analyzed in order to develop "Universal" 28S probes. All sequences were subjected to extensive manipulation to identify optimal relative

alignments in order to identify regions of similarity for use as "Universal" probes. The selected probe sequences had to meet several important criteria besides the condition of being present in 28S genes from most fungal species. Each probe sequence required an appropriate thermal profile, secondary structure and utility in a DNA amplification reaction. These probes were optimized to work for PCR amplification in pure cultures of fungus, as well as in the presence of DNA from multiple sources as in the case of clinical specimens. The probes were also designed to facilitate direct sequencing of the amplified DNA. Our analysis led to the discovery of the oligonucleotide probes listed in (SEQ ID NO:1) and (SEQ ID NO:2). (For their location, see Figure 1.) The successful identification of these two probes ((SEQ ID NO:1) and (SEQ ID NO:2)) completed our first objective to develop nucleic acid probes that would hybridize to, and detect 28S rRNA and rDNA from a majority of fungi (Figure 1 and Table 1). As shown later in this application, the novel sequence information generated by the use of our "Universal" probes allowed us to develop species-specific probes ((SEQ ID NO:3) to (SEQ ID NO:23)) capable of identifying 19 different disease-causing fungi.

Table 1:

Presence of hybridization sites for probes SEQ ID NO: 1 and SEQ ID NO: 2 in 28S nucleic acid sequences.

	SEQ ID NO: 1	SEQ ID NO: 2
<i>Acremonium</i> sp.	+	+
<i>Aspergillus clavatus</i>	+	+
<i>Aspergillus flavus</i>	+	+
<i>Aspergillus fumigatus</i>	+	+
<i>Aspergillus glaucus</i>	+	+
<i>Aspergillus nidulans</i>	+	+
<i>Aspergillus niger</i>	+	+
<i>Aspergillus ochraceus</i>	+	+

<i>Aspergillus terreus</i>	+	+
<i>Aspergillus unguis</i>	+	+
<i>Aspergillus ustus</i>	+	+
<i>Beauveria</i> sp.	+	+
<i>Bipolaris</i> sp.	+	+
<i>Blastomyces dermatitidis</i>	+	+
<i>Blastoschizomyces</i> sp.	+	+
<i>Candida albicans</i>	+	+
<i>Candida glabrata</i>	+	+
<i>Candida guilliermondii</i>	+	+
<i>Candida kefyr</i>	+	+
<i>Candida krusei</i>	+	+
<i>Candida lusitanae</i>	+	+
<i>Candida parapsilosis</i>	+	+
<i>Candida tropicalis</i>	+	+
<i>Chrysosporium</i> sp.	+	+
<i>Cladosporium</i> sp.	+	+
<i>Coccidioides immitis</i>	+	+
<i>Cryptococcus laurentii</i>	+	+
<i>Cryptococcus neoformans</i> serotype A	+	+
<i>Cryptococcus neoformans</i> var. <i>gattii</i> serotype B	+	+
<i>Cryptococcus terreus</i>	+	+
<i>Curvularia</i> sp.	+	+
<i>Filobasidiella (Cryptococcus) neoformans</i> var <i>bacillispora</i> serotype C	+	+
<i>Filobasidiella (Cryptococcus) neoformans</i> var <i>neoformans</i> serotype D	+	+
<i>Filobasidium capsuligenum</i>	+	+
<i>Filobasidium uniguttulatum</i>	+	+
<i>Fusarium</i> sp.	+	+
<i>Geotrichum</i> sp.	+	+
<i>Histoplasma capsulatum</i>	+	+
<i>Malbranchea</i> sp.	+	+
<i>Mucor</i> sp.	+	+
<i>Paecilomyces</i> sp.	+	+
<i>Penicillium</i> sp.	+	+
<i>Pseudallescheria boydii</i>	+	+
<i>Rhizopus</i> sp.	+	+



<i>Saccharomyces cerevisiae</i>	+	+
<i>Scopulariopsis brevicaulis</i>	+	+
<i>Scopulariopsis brumptii</i>	+	+
<i>Sporothrix schenckii</i>	+	+
<i>Trichosporon beigelia</i>	+	+
Human	-	+

- Probes SEQ ID NO: 1 and SEQ ID NO: 2 were used to successfully amplify (Table 2) and sequence DNA (Figure 2) spanning this variable region from the following 49 organisms: *Acremonium* sp., *Aspergillus clavatus*, *Aspergillus flavus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*, *Aspergillus ustus*, *Beauveria* sp., *Bipolaris* sp., *Blastomyces dermatitidis*, *Blastoschizomyces* sp., *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* serotype A, *Cryptococcus neoformans* var. *gattii* serotype B, *Cryptococcus terreus*, *Cryptococcus laurentii*, *Curvularia* sp., *Filobasidiella* (*Cryptococcus*) *neoformans* var. *bacillispora* serotype C, *Filobasidiella* (*Cryptococcus*) *neoformans* var. *neoformans* serotype D, *Filobasidium capsuligenum*, *Filobasidium uniguttulatum*, *Fusarium* sp., *Geotrichum* sp., *Histoplasma capsulatum*, *Malbranchea* sp., *Mucor* sp., *Paecilomyces* sp., *Penicillium* sp., *Pseudallescheria boydii*, *Rhizopus* sp., *Saccharomyces cerevisiae*, *Scopulariopsis brevicaulis*, *Scopulariopsis brumptii*, *Sporothrix schenckii* and *Trichosporon beigelia*. This list contains all 4 serotypes (A, B, C and D) of *Cryptococcus neoformans*. This sequence information generated by the use of probes SEQ ID NO: 1 and SEQ ID NO: 2 expanded the size of our database consisting of fungal 28S sequences. All amplified DNA was sequenced across both strands from a minimum of two different isolates of each organism to ensure accuracy of the data generated.

Table 2:

Polymerase chain reaction amplification of 28S rDNA with probes SEQ ID NO: 1  
and

5 SEQ ID NO: 2.

	PCR with SEQ ID NO: 1 & NO: 2
<i>Acremonium</i> sp.	+
<i>Aspergillus clavatus</i>	+
<i>Aspergillus flavus</i>	+
<i>Aspergillus fumigatus</i>	+
<i>Aspergillus glaucus</i>	+
<i>Aspergillus nidulans</i>	+
<i>Aspergillus niger</i>	+
<i>Aspergillus ochraceus</i>	+
<i>Aspergillus terreus</i>	+
<i>Aspergillus unguis</i>	+
<i>Aspergillus ustus</i>	+
<i>Beauveria</i> sp.	+
<i>Bipolaris</i> sp.	+
<i>Blastomyces dermatitidis</i>	+
<i>Blastoschizomyces</i> sp.	+
<i>Candida albicans</i>	+
<i>Candida glabrata</i>	+
<i>Candida guilliermondii</i>	+
<i>Candida kefyr</i>	+
<i>Candida krusei</i>	+
<i>Candida lusitanae</i>	+
<i>Candida parapsilosis</i>	+
<i>Candida tropicalis</i>	+
<i>Chrysosporium</i> sp.	+
<i>Cladosporium</i> sp.	+
<i>Coccidioides immitis</i>	+
<i>Cryptococcus laurentii</i>	+
<i>Cryptococcus neoformans</i> serotype A	+
<i>Cryptococcus neoformans</i> var. <i>gattii</i> serotype B	+
<i>Cryptococcus terreus</i>	+

<i>Curvularia</i> sp.	+
<i>Filobasidiella</i> ( <i>Cryptococcus</i> ) <i>neoformans</i> var <i>bacillispora</i> serotype C	+
<i>Filobasidiella</i> ( <i>Cryptococcus</i> ) <i>neoformans</i> var <i>neoformans</i> serotype D	+
<i>Filobasidium capsuligenum</i>	+
<i>Filobasidium uniguttulatum</i>	+
<i>Fusarium</i> sp.	+
<i>Geotrichum</i> sp.	+
<i>Histoplasma capsulatum</i>	+
<i>Malbranchea</i> sp.	+
<i>Mucor</i> sp.	+
<i>Paecilomyces</i> sp.	+
<i>Penicillium</i> sp.	+
<i>Pseudallescheria boydii</i>	+
<i>Rhizopus</i> sp.	+
<i>Saccharomyces cerevisiae</i>	+
<i>Scopulariopsis brevicaulis</i>	+
<i>Scopulariopsis brumptii</i>	+
<i>Sporothrix schenckii</i>	+
<i>Trichosporon beigelii</i>	+
Human	-

This list of fungi sequenced by us represents organisms responsible for most cases of subcutaneous and deep mycotic infections in humans and also includes saprophytes (non-pathogenic fungi) commonly encountered in clinical isolates. Since the two probes (SEQ ID NO: 1 and SEQ ID NO: 2) hybridize to 28S rDNA from all the fungi listed above, they are capable of diagnosing the presence of a majority of fungi that are likely to be present in a clinical specimen. They are believed to be primers for universally detecting fungi.

10

Probes listed in SEQ ID NO: 1 and SEQ ID NO: 2 were also checked for their potential ability to hybridize to, and amplify (in a polymerase chain reaction) 23S sequences from bacteria by searching for hybridization sites among the 539

bacterial 23S genes listed in GenBank. Bacterial 23S rDNAs do not have suitable hybridization sites for SEQ ID NO: 1 and SEQ ID NO: 2 and these two probes should not be able to amplify bacterial DNA under stringent conditions.

- 5           Our second objective was to develop species specific probes, which under hybridization conditions, would detect *Aspergillus fumigatus*, *Blastomyces dermatitidis*, *Candida albicans*, *Coccidioides immitis*, *Cryptococcus neoformans*, *Histoplasma capsulatum*, *Aspergillus flavus*, *Aspergillus glaucus*, *Aspergillus niger*, *Aspergillus terreus*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*,  
10 *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Pseudallescheria boydii*, and *Sporothrix schenckii*. We used our database of fungal 28S nucleic acid sequences to create a multiple sequence alignment of all the organisms that we had sequenced. Every individual sequence was subjected to intensive comparison with all other sequences in our database in order to discover  
15 unique regions of sequence that would be present only in the fungus of interest, and would be absent in all other fungi. When unique stretches of sequence were identified, these were further analyzed for thermal profile and secondary structure. Each probe constructed by us will, under conditions of hybridization, specifically hybridize to and detect, nucleic acid sequence from the unique region of only one  
20 specific target fungus. Those versed in the art will recognize that specification of a single-stranded DNA sequence implies the utility of the complementary DNA sequence, as well as the two equivalent RNA sequences. Furthermore, sequences incorporating modification of any of the moieties comprising the nucleic acid (i.e., the base, the sugar or the backbone) are functional equivalents of the sequence. It  
25 should also be recognized that these additional sequences can potentially serve as probes or primers. Finally, those versed in the art recognize that comparisons of extensive DNA sequences provides enough variability and uniqueness to speciate organisms (Figure 2).

The nucleic acid sequences for these species specific synthetic probes are listed in SEQ ID NO: 3 to SEQ ID NO: 23. There are two probes specific for *Cryptococcus neoformans*, two probes specific for *Sporothrix schenckii*, and one probe each for *Aspergillus fumigatus*, *Blastomyces dermatitidis*, *Candida albicans*,  
 5 *Coccidioides immitis*, *Histoplasma capsulatum*, *Aspergillus flavus*, *Aspergillus glaucus*, *Aspergillus niger*, *Aspergillus terreus*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis* and *Pseudallescheria boydii* 28S rRNA and rDNA. (See Tables 3 - 6 and further discussion below.)

10

All species specific probes developed by us are novel and to the best of our knowledge have not been reported in the literature. While all 28S genes sequenced by us had several regions that were different among the various species analyzed, the regions that would function best as species specific probes under conditions of  
 15 hybridization were not obvious. Extensive analysis of each 28S sequence yielded several potential probe sites. These were studied in detail to enable the selection of optimal unique sites for each probe, based on the need to obtain optimal hybridization characteristics under the test conditions. The highly specific hybridization characteristics of all probe sequences developed by us were then  
 20 validated by experimental results. The prior existence in GenBank of sequences for *Candida albicans* and only one serotype of *Cryptococcus neoformans* 28S genes was in itself not sufficient to enable even an individual versed in this field to develop specific probes for either of these two organisms. We had to obtain novel 28S sequence from *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*,  
 25 *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Cryptococcus neoformans* serotype A, *Cryptococcus neoformans* var. *gattii* serotype B, *Cryptococcus terreus*, *Cryptococcus laurentii*, *Filobasidiella* (*Cryptococcus*) *neoformans* var *bacillispora* serotype C, *Filobasidiella* (*Cryptococcus*) *neoformans* var *neoformans* serotype D, *Filobasidium*

*capsuligenum* and *Filobasidium uniguttulatum* before we were able to identify potential regions for the development of species specific probes for these two fungal organisms that would not cross react with the others listed above.

5           Our modification of the Chomczynski technique (see Example 2, below) allows us to obtain DNA from any clinical specimen, irrespective of source (see Table 8 for a variety of clinical specimens tested), within a 3 hour period. The PCR amplification and subsequent probing can be accomplished with ease within a 24  
10           hour period. The final identification is therefore possible in a day as opposed to several days or weeks required by traditional methods. This speed and sensitivity of  
15           diagnosis can make a difference between life and death in debilitated patients battling fungal diseases of undetermined cause. Rapid diagnosis will allow physicians to immediately direct their therapy towards curing the identified causative fungus, rather than wait for days or weeks while the patient succumbs to  
15           an unknown fungus.

          Our probes have the ability to pick out the correct target organism even in a mixed fungal infection because of their high level of specificity. The methods of Hopfer et. al. and Maiwald et. al., do not allow identification of individual species  
20           in a mixed fungal infection because restriction fragment length polymorphism results are nearly impossible to interpret when multiple organisms contribute to the restriction fragments. Their method can therefore only be used on a pure culture, and this also does not save any diagnostic time, because the fungus first has to be  
25           grown in culture.

          The probes developed by us allow rapid species identification of a large number of pathogenic fungi by using multiple probes against only one PCR amplified fragment of DNA. Coupled with our modified DNA extraction technique and our ability to accurately diagnose in the case of mixed organisms, this strategy

can provide the greatest amount of diagnostic information in the shortest amount of time. This diagnostic strategy is also amenable to automation, which can result in even greater savings in time, money and effort.

5           The sequences and the complement of the sequences claimed in this disclosure, along with any modifications to these sequences, may potentially be utilized in assays for the identification of fungi based on several existing methodologies, as well as future improvements and alterations of this technology. These techniques include, but are not limited to, assays based on hybridization, 10 ligation, polymerization, depolymerization, sequencing, chemical degradation, enzymatic digestion, electrophoresis, chromatography and amplification. Furthermore, all such variations ultimately are based in some selection or amplification process, some ligand or some nucleic acid moiety that recognizes or utilizes the sequences (SEQ ID NO: 1) to (SEQ ID NO:23) claimed in this 15 application. Such variations include but are not limited to use of a variety of linear or exponential target amplification schemes, such as, any of the myriad forms of PCR, the ligase chain reaction, Q-beta repliase, etc.; direct detection of species-specific nucleic acid purified or extracted from pure fungal culture using a probe selected from the group (SEQ ID NO: 3) to (SEQ ID NO: 23); use of the 20 complementary DNA forms of (SEQ ID NO:1) to (SEQ ID NO:23); use of the RNA forms of these sequences and their complements; and use of derivatives of these DNA or RNA sequences by the addition of one or more reporter moieties from a variety of labels including nucleic acid sequences, proteins, signal generating ligands such as acridinium esters, and/or paramagnetic particles. These techniques 25 may be utilized with DNA, RNA or modified derivatives used as either the target or the detection molecule.

In addition to the 23 sequences SEQ ID NO: 1 to SEQ ID NO: 23, we also describe an additional 51 sequences SEQ ID NO: 24 to SEQ ID NO: 74. These 51

sequences are inclusive of SEQ ID NO: 3 to SEQ ID NO: 23 and are shown as a multiple sequence alignment (Figure 2) with coordinate 1 corresponding to base # 431 of a reference *S. cerevisiae* 28S rRNA gene. (The numbers are comparable to the primary sequence of *S. cerevisiae* 28S rRNA gene. Genbank accession number: 5 J01355). These sequences were obtained by amplifying and sequencing 28S rDNA from various fungi with primers SEQ ID NO: 1 and SEQ ID NO: 2. (SEQ ID NO: 1 corresponds to coordinates 403-422 and the SEQ ID NO: 2 corresponds to coordinates 645-662 of the reference *S. cerevisiae* gene).

10       An analysis of these aligned sequences enabled us to develop the species specific probes SEQ ID NO: 3 to SEQ ID NO: 23, and sites for these probes are shown underlined. These 51 aligned sequences contain sufficient variability, to enable a person versed in this art, to develop additional species specific hybridization probes in the 10-50 nucleotide length. Similarly, longer species 15 specific hybridization probes encompassing the entire 200+ nucleotide length can also be envisioned. Species identification may also be accomplished by direct DNA sequence determination of any DNA amplified with primers SEQ ID NO: 1 and SEQ ID NO: 2. If the derived sequence matches approximately 98% or more of any sequence in SEQ ID NO: 24 to SEQ ID NO: 74, then the identity of the organism 20 can be ascertained. Additionally, we recognize that parts of SEQ ID NO: 24 to SEQ ID NO: 74 may be specific for groups of fungi arranged phylogenetically at the level of genus or higher. SEQ ID NO: 24 to SEQ ID NO: 74, their complements, along with any modification to these sequences may also potentially be utilized in assays for the identification of fungi based on existing methodologies 25 and future technologies as noted above for SEQ ID NO: 1 to SEQ ID NO: 23.



## Legend to figure 2:

The multiple sequence alignment shows the sequence of 28S ribosomal RNA genes amplified with primers SEQ ID NO: 1 and SEQ ID NO: 2. 21 species specific probes (SEQ ID NO: 3 to SEQ ID NO: 23) are shown underlined. Minor sequence variation among two isolate of the same organism are represented by the appropriate code (see key below). Major differences among *Rhizopus* species are depicted by including 3 separate *Rhizopus* sequences in the alignment. (The organisms in this figure are listed according to their sequence relatedness.)

10

## Key to symbols:

(.) gap in sequence to facilitate alignment

(R) A or G

15 (W) A or T

(Y) T or C

(M) A or C

(K) T or G

(S) G or C

20 (B) T,G or C

Acremo *Acremonium* species

A\_clav *Aspergillus clavatus*

A\_flav *Aspergillus flavus*

25 A\_fumi *Aspergillus fumigatus*

A\_glau *Aspergillus glaucus*

A\_nidu *Aspergillus nidulans*

A\_nige *Aspergillus niger*

A\_ochr *Aspergillus ochraceus*

	A_terr	<i>Aspergillus terreus</i>
	A_ungu	<i>Aspergillus unguis</i>
	A_ustu	<i>Aspergillus ustus</i>
	Beauve	<i>Beauveria</i> species
5	Bipola	<i>Bipolaris</i> species
	Blasch	<i>Blastoschizomyces</i> species
	B_derm	<i>Blastomyces dermatitidis</i>
	Chryso	<i>Chrysosporium</i> species
	Clados	<i>Cladosporium</i> species
10	Curvul	<i>Curvularia</i> species
	C_albi	<i>Candida albicans</i>
	C_glab	<i>Candida glabrata</i>
	C_guil	<i>Candida guilliermondii</i>
	C_immi	<i>Coccidioides immitis</i>
15	C_kefy	<i>Candida kefyr</i>
	C_krus	<i>Candida krusei</i>
	C_laur	<i>Cryptococcus laurentii</i>
	C_lusi	<i>Candida lusitanae</i>
	C_neob	<i>Cryptococcus neoformans</i> var <i>gattii</i> serotype B
20	C_neof	<i>Cryptococcus neoformans</i> serotype A
	C_para	<i>Candida parapsilosis</i>
	C_terr	<i>Cryptococcus terreus</i>
	C_trop	<i>Candida tropicalis</i>
	Fusari	<i>Fusarium</i> species
25	F_caps	<i>Filobasidium capsuligenum</i>
	F_neoc	<i>Filobasidiella (Cryptococcus) neoformans</i> var <i>bacillispora</i> serotype C
	F_neod	<i>Filobasidiella (Cryptococcus) neoformans</i> var <i>neoformans</i> serotype D

	F_unig	<i>Filobasidium uniguttulatum</i>
	Geotri	<i>Geotrichum</i> species
	H_caps	<i>Histoplasma capsulatum</i>
	Malbra	<i>Malbranchea</i> species
5	Mucor_	<i>Mucor</i> species
	Paecil	<i>Paecilomyces</i> species
	Penici	<i>Penicillium</i> species
	P_boyd	<i>Pseudallescheria boydii</i>
	Rhizo1	<i>Rhizopus</i> species isolate #1
10	Rhizo2	<i>Rhizopus</i> species isolate #2
	Rhizo3	<i>Rhizopus</i> species isolate #3
	Sporot	<i>Sporothrix schenkii</i>
	S_brev	<i>Scopulariopsis brevicaulis</i>
	S_brum	<i>Scopulariopsis brumpti</i>
15	S_cere	<i>Saccharomyces cerevisiae</i>
	T_beig	<i>Trichosporon beigelii</i>

Further variations of the invention that utilize any of the named sequences  
 20 will be apparent to those with ordinary skill in the art. The following examples  
 illustrate various aspects of the invention but are not intended to limit its usefulness.

EXAMPLE 1. Testing probes SEQ ID NO: 3 to SEQ ID NO: 23 for hybridization  
 specificity.

25

Probes listed in SEQ ID NO: 3 to SEQ ID NO: 23 were tested for specificity  
 against their target organisms. Probe SEQ ID NO: 5 for *Candida albicans* was the  
 first one tested against a panel of fungi taken from the Mayo Clinic collection. 28S  
 rDNA from *Acremonium* sp., *Aspergillus clavatus*, *Aspergillus flavus*, *Aspergillus*

*fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*, *Aspergillus ustus*, *Aspergillus* sp., *Beauveria* sp., *Bipolaris* sp., *Blastomyces dermatitidis*, *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*,  
5 *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* serotype A, *Curvularia* sp., *Fusarium* sp., *Geotrichum* sp., *Histoplasma capsulatum*, *Mucor* sp., *Penicillium* sp., *Pseudallescheria boydii*, *Rhizopus* sp., *Saccharomyces cerevisiae*, *Scopulariopsis brevicaulis*, *Sporothrix schenckii* and *Trichosporon beigeli* was  
10 amplified in a polymerase chain reaction using oligonucleotide probes SEQ ID NO: 1 and SEQ ID NO: 2. All PCR amplifications were carried out as hot-start reactions in a 50 ul reaction volume using Perkin-Elmer (Norwalk, CT) 0.5 ml thin-wall polypropylene tubes and a Perkin-Elmer thermal cycler. Reagents added to the tube initially were 2.5 ul of 10X PCR buffer (100 mM tris pH 8.3, 500 mM KCl, 15  
15 mM MgCl<sub>2</sub>), 5.0 ul of 50% glycerol/1 mM cresol red, 8.0 ul of dNTP mix (1.25 mM each of dATP, dGTP, dTTP and dCTP), 11 picomoles of each nucleic acid primer and sterile water to make up a volume of 25 ul. A wax bead (Ampliwax Gem=100, Perkin-Elmer) was added and the tubes heated to 77°C for 30 seconds and cooled to room temperature to form a wax barrier. 2.5 ul of 10X PCR buffer,  
20 5.0 ul of 50% glycerol/1 mM cresol red, 0.2 ul Taq polymerase (AmpliTaq 5U/ul, Perkin-Elmer) and 15.3 ul of sterile water was added to the tube along with 2.0 ul of DNA from the fungal whole cell boiled lysate described above. 50 cycles of thermal cycling was carried out at 94°C - 30 sec, 50°C - 1 min, 72°C - 2 min. Five microliters of polymerase chain reaction mix from each sample was run on a 5%  
25 polyacrylamide gel to visually confirm the successful amplification of 28S rDNA from each fungus listed above. 40 ul of the remaining amplified 28S rDNA was denatured in 1 N NaOH, and half of this denatured rDNA was slot blotted on to a positively charged polysulphone based membrane equilibrated in 0.5 N NaOH. The membrane was air dried for 15 minutes and baked in a vacuum oven at 80°C for 30

minutes. Amplified rDNA from each species was now bound and immobilized at a separate spot on the membrane. The free binding sites on the membrane were blocked by incubating the membrane for 3 hours at 40°C in hybridization buffer (100 ml of hybridization buffer was made using 1g non-fat milk powder, 6g NaH<sub>2</sub>PO<sub>4</sub>, 7g SDS, 200 ul 0.5M EDTA and adjusted to pH 7.2 with NaOH). The specific probe for *Candida albicans* (SEQ ID NO: 5) was end-labeled with radioactive phosphorus using <sup>32</sup>P ATP and T4 polynucleotide kinase. 50 picomoles of this probe was added to 70 milliliters of hybridization buffer and the membrane was probed at 40°C overnight. The membrane was washed in hybridization buffer at 40°C for 15 minutes followed by a wash in 2X SSC at 40°C for 15 minutes. The membrane was then exposed on x-ray film for at least 1 hour. The oligonucleotide probe SEQ ID NO: 5 only hybridized to amplified 28S rDNA from *Candida albicans* (see Table 3) Under these hybridization conditions, probe SEQ ID NO: 5 is extremely specific for *Candida albicans*. The sequence of oligonucleotide probe SEQ ID NO: 5 differs from the sequences of other species of *Candida* by as few as 1 or 2 bases, but these mismatches are sufficient to prevent stable hybrids from forming with the other *Candida* species.

Probes SEQ ID NO: 3 to SEQ ID NO: 23 were tested for specificity, as described above for the *Candida albicans* probe SEQ ID NO: 5, against the same panel of fungi listed in the preceding paragraph. The positively charged polysulphone based membrane probed with *Candida albicans* probe SEQ ID NO: 5 was washed in 0.5 N NaOH at 40°C for 10 minutes to remove all bound *Candida albicans* probe. The membrane was sequentially probed with all probes listed in SEQ ID NO: 3 to SEQ ID NO: 23. For each subsequently tested probe, the membrane was blocked for at least 30 minutes, probe hybridization was carried out at 40-42°C for at least 3 hours, and post-hybridization washes were done in 2X SSC for 20 minutes. The membrane was stripped between probings by washing in 0.5 to 1.0 N NaOH at 40-42°C. Results are listed in Tables 3 to 6.

As shown in Tables 3 to 6, each probe listed in SEQ ID NO: 3 to SEQ ID NO: 23 specifically hybridizes to only one target fungal 28S nucleic acid sequence. This specificity is essential for identifying a given species of fungus in clinical specimens containing mixed fungal organisms with a high level of reliability. The

5 39 organisms listed in these Tables represent a majority of organisms that are commonly isolated from clinical samples. While we have developed 21 species specific probes (SEQ ID NO: 3 to SEQ ID NO: 23) that identify a total of 19 individual organisms, the additional organisms listed in the test panel were used to ensure that our probes did not have any cross-reactivity with other fungi likely to be

10 present in a clinical specimen. The ability to accurately and reliably diagnose, and identify to a species level, this large a number of pathogens is unmatched by any other report. The fact that we can achieve this by probing DNA amplified by a single pair of "Universal" probes (SEQ ID NO: 1 and SEQ ID NO: 2) is highly advantageous as it saves time, money and effort by providing the ability to test a

15 single amplified target with 21 different probes (SEQ ID NO: 3 to SEQ ID NO: 23).

A GenBank search was carried out with all probes listed in SEQ ID NO: 3 to SEQ ID NO: 23 in order to determine whether similar gene sequences were present in the database. 28S sequences for *Candida albicans* and one serotype of

20 *Cryptococcus neoformans* are already present in GenBank, and as expected, the probes for *Candida albicans* and *Cryptococcus neoformans* correctly identified the 28S sequences from these two organisms. Ten other probes also matched DNA sequences from a variety of genes not related to the 28S gene (Table 7). This was expected because short stretches of sequence identity can often be found for any

25 query sequence in unrelated genes from the same or a different organism. This observation is known to those versed in this art. In all cases, sequences that matched a probe sequence were not located within the 28S rRNA genes. Our probes are used to analyze 28S DNA that has been previously amplified in a polymerase chain reaction with our probes SEQ ID NO: 1 and SEQ ID NO: 2. Under stringent

conditions, these two probes only amplify DNA from fungal 28S rRNA genes. Therefore no amplified DNA from the non-28S genes listed in Table 7 will be available for the hybridization of probes SEQ ID NO: 3 to SEQ ID NO: 23. The presence of related sequences in non-28S, unamplified genes will not be detected  
5 and will, thus, not have any effect on the sensitivity or the specificity of our detection and identification strategy.

Table 3:

Detection of species specific 28S sequence with probes SEQ ID NO: 3 to SEQ ID NO: 8

FUNGUS	SEQ ID: 3	SEQ ID: 4	SEQ ID: 5	SEQ ID: 6	SEQ ID: 7	SEQ ID: 8
<i>Acremonium</i> sp.	-	-	-	-	-	-
<i>Aspergillus clavatus</i>	-	-	-	-	-	-
<i>Aspergillus flavus</i>	-	-	-	-	-	-
<i>Aspergillus fumigatus</i>	+	-	-	-	-	-
<i>Aspergillus glaucus</i>	-	-	-	-	-	-
<i>Aspergillus nidulans</i>	-	-	-	-	-	-
<i>Aspergillus niger</i>	-	-	-	-	-	-
<i>Aspergillus ochraceus</i>	-	-	-	-	-	-
<i>Aspergillus terreus</i>	-	-	-	-	-	-
<i>Aspergillus unguis</i>	-	-	-	-	-	-
<i>Aspergillus ustus</i>	-	-	-	-	-	-
<i>Aspergillus</i> sp.	-	-	-	-	-	-
<i>Beauveria</i> sp.	-	-	-	-	-	-
<i>Bipolaris</i> sp.	-	-	-	-	-	-
<i>Blastomyces dermatitidis</i>	-	+	-	-	-	-
<i>Candida albicans</i>	-	-	+	-	-	-
<i>Candida glabrata</i>	-	-	-	-	-	-
<i>Candida guilliermondii</i>	-	-	-	-	-	-
<i>Candida kefyr</i>	-	-	-	-	-	-
<i>Candida krusei</i>	-	-	-	-	-	-
<i>Candida lusitanae</i>	-	-	-	-	-	-
<i>Candida parapsilosis</i>	-	-	-	-	-	-
<i>Candida tropicalis</i>	-	-	-	-	-	-
<i>Chrysosporium</i> sp.	-	-	-	-	-	-
<i>Cladosporium</i> sp.	-	-	-	-	-	-
<i>Coccidioides immitis</i>	-	-	-	+	-	-
<i>Cryptococcus neoformans</i>	-	-	-	-	+	+
<i>Curvularia</i> sp.	-	-	-	-	-	-
<i>Fusarium</i> sp.	-	-	-	-	-	-
<i>Geotrichum</i> sp.	-	-	-	-	-	-
<i>Histoplasma capsulatum</i>	-	-	-	-	-	-
<i>Mucor</i> sp.	-	-	-	-	-	-
<i>Penicillium</i> sp.	-	-	-	-	-	-
<i>Pseudallescheria boydii</i>	-	-	-	-	-	-
<i>Rhizopus</i> sp.	-	-	-	-	-	-
<i>Saccharomyces cerevisiae</i>	-	-	-	-	-	-
<i>Scopulariopsis brevicaulis</i>	-	-	-	-	-	-
<i>Sporothrix schenckii</i>	-	-	-	-	-	-
<i>Trichosporon beigeli</i>	-	-	-	-	-	-



+	Positive
-	Negative after 20 minute wash in 2X SSC

Table 4:

Detection of species specific 28S sequence with probes SEQ ID NO: 9 to SEQ ID NO: 14

5

FUNGUS	SEQ ID: 9	SEQ ID: 10	SEQ ID: 11	SEQ ID: 12	SEQ ID: 13	SEQ ID: 14
<i>Acremonium</i> sp.	-	-	-	-	-	-
<i>Aspergillus clavatus</i>	-	-	-	-	-	-
<i>Aspergillus flavus</i>	-	-	-	-	-	-
<i>Aspergillus fumigatus</i>	-	-	-	-	-	-
<i>Aspergillus glaucus</i>	-	+	-	-	-	-
<i>Aspergillus nidulans</i>	-	-	-	-	-	-
<i>Aspergillus niger</i>	-	-	+	-	-	-
<i>Aspergillus ochraceus</i>	-	-	-	-	-	-
<i>Aspergillus terreus</i>	-	-	-	+	-	-
<i>Aspergillus unguis</i>	-	-	-	-	-	-
<i>Aspergillus ustus</i>	-	-	-	-	-	-
<i>Aspergillus</i> sp.	-	-	-	-	-	-
<i>Beauveria</i> sp.	-	-	-	-	-	-
<i>Bipolaris</i> sp.	-	-	-	-	-	-
<i>Blastomyces dermatitidis</i>	-	-	-	-	-	-
<i>Candida albicans</i>	-	-	-	-	-	-
<i>Candida glabrata</i>	-	-	-	-	+	-
<i>Candida guilliermondii</i>	-	-	-	-	-	+
<i>Candida kefyr</i>	-	-	-	-	-	-
<i>Candida krusei</i>	-	-	-	-	-	-
<i>Candida lusitanae</i>	-	-	-	-	-	-
<i>Candida parapsilosis</i>	-	-	-	-	-	-
<i>Candida tropicalis</i>	-	-	-	-	-	-
<i>Chrysosporium</i> sp.	-	-	-	-	-	-
<i>Cladosporium</i> sp.	-	-	-	-	-	-
<i>Coccidioides immitis</i>	-	-	-	-	-	-
<i>Cryptococcus neoformans</i>	-	-	-	-	-	-
<i>Curvularia</i> sp.	-	-	-	-	-	-
<i>Fusarium</i> sp.	-	-	-	-	-	-
<i>Geotrichum</i> sp.	-	-	-	-	-	-
<i>Histoplasma capsulatum</i>	+	-	-	-	-	-
<i>Mucor</i> sp.	-	-	-	-	-	-
<i>Penicillium</i> sp.	-	-	-	-	-	-
<i>Pseudallescheria boydii</i>	-	-	-	-	-	-
<i>Rhizopus</i> sp.	-	-	-	-	-	-
<i>Saccharomyces cerevisiae</i>	-	-	-	-	-	-
<i>Scopulariopsis brevicaulis</i>	-	-	-	-	-	-
<i>Sporothrix schenckii</i>	-	-	-	-	-	-
<i>Trichosporon beigeli</i>	-	-	-	-	-	-

+	Positive
-	Negative after 20 minute wash in 2X SSC

## 5 Table 5:

Detection of species specific 28S sequence with probes SEQ ID NO: 15 to SEQ ID NO: 20

FUNGUS	SEQ ID: 15	SEQ ID: 16	SEQ ID: 17	SEQ ID: 18	SEQ ID: 19	SEQ ID: 20
<i>Acremonium</i> sp.	-	-	-	-	-	-
<i>Aspergillus clavatus</i>	-	-	-	-	-	-
<i>Aspergillus flavus</i>	-	-	-	-	-	-
<i>Aspergillus fumigatus</i>	-	-	-	-	-	-
<i>Aspergillus glaucus</i>	-	-	-	-	-	-
<i>Aspergillus nidulans</i>	-	-	-	-	-	-
<i>Aspergillus niger</i>	-	-	-	-	-	-
<i>Aspergillus ochraceus</i>	-	-	-	-	-	-
<i>Aspergillus terreus</i>	-	-	-	-	-	-
<i>Aspergillus unguis</i>	-	-	-	-	-	-
<i>Aspergillus ustus</i>	-	-	-	-	-	-
<i>Aspergillus</i> sp.	-	-	-	-	-	-
<i>Beauveria</i> sp.	-	-	-	-	-	-
<i>Bipolaris</i> sp.	-	-	-	-	-	-
<i>Blastomyces dermatitidis</i>	-	-	-	-	-	-
<i>Candida albicans</i>	-	-	-	-	-	-
<i>Candida glabrata</i>	-	-	-	-	-	-
<i>Candida guilliermondii</i>	-	-	-	-	-	-
<i>Candida kefyr</i>	+	-	-	-	-	-
<i>Candida krusei</i>	-	+	-	-	-	-
<i>Candida lusitanae</i>	-	-	+	-	-	-
<i>Candida parapsilosis</i>	-	-	-	+	-	-
<i>Candida tropicalis</i>	-	-	-	-	+	-
<i>Chrysosporium</i> sp.	-	-	-	-	-	-
<i>Cladosporium</i> sp.	-	-	-	-	-	-
<i>Coccidioides immitis</i>	-	-	-	-	-	-
<i>Cryptococcus neoformans</i>	-	-	-	-	-	-
<i>Curvularia</i> sp.	-	-	-	-	-	-
<i>Fusarium</i> sp.	-	-	-	-	-	-
<i>Geotrichum</i> sp.	-	-	-	-	-	-
<i>Histoplasma capsulatum</i>	-	-	-	-	-	-
<i>Mucor</i> sp.	-	-	-	-	-	-
<i>Penicillium</i> sp.	-	-	-	-	-	-
<i>Pseudallescheria boydii</i>	-	-	-	-	-	+
<i>Rhizopus</i> sp.	-	-	-	-	-	-
<i>Saccharomyces cerevisiae</i>	-	-	-	-	-	-

<i>Scopulariopsis brevicaulis</i>	-	-	-	-	-	-
<i>Sporothrix schenckii</i>	-	-	-	-	-	-
<i>Trichosporon beigeli</i>	-	-	-	-	-	-

+	Positive
-	Negative after 20 minute wash in 2X SSC

## 5 Table 6:

Detection of species specific 28S sequence with probes SEQ ID NO: 21 to SEQ ID NO: 23

FUNGUS	SEQ ID: 21	SEQ ID: 22	SEQ ID: 23
<i>Acremonium</i> sp.	-	-	-
<i>Aspergillus clavatus</i>	-	-	-
<i>Aspergillus flavus</i>	+	-	-
<i>Aspergillus fumigatus</i>	-	-	-
<i>Aspergillus glaucus</i>	-	-	-
<i>Aspergillus nidulans</i>	-	-	-
<i>Aspergillus niger</i>	-	-	-
<i>Aspergillus ochraceus</i>	-	-	-
<i>Aspergillus terreus</i>	-	-	-
<i>Aspergillus unguis</i>	-	-	-
<i>Aspergillus ustus</i>	-	-	-
<i>Aspergillus</i> sp.	-	-	-
<i>Beauveria</i> sp.	-	-	-
<i>Bipolaris</i> sp.	-	-	-
<i>Blastomyces dermatitidis</i>	-	-	-
<i>Candida albicans</i>	-	-	-
<i>Candida glabrata</i>	-	-	-
<i>Candida guilliermondii</i>	-	-	-
<i>Candida kefyr</i>	-	-	-
<i>Candida krusei</i>	-	-	-
<i>Candida lusitanae</i>	-	-	-
<i>Candida parapsilosis</i>	-	-	-
<i>Candida tropicalis</i>	-	-	-
<i>Chrysosporium</i> sp.	-	-	-
<i>Cladosporium</i> sp.	-	-	-
<i>Coccidioides immitis</i>	-	-	-
<i>Cryptococcus neoformans</i>	-	-	-
<i>Curvularia</i> sp.	-	-	-
<i>Fusarium</i> sp.	-	-	-
<i>Geotrichum</i> sp.	-	-	-
<i>Histoplasma capsulatum</i>	-	-	-
<i>Mucor</i> sp.	-	-	-
<i>Penicillium</i> sp.	-	-	-

<i>Pseudallescheria boydii</i>	-	-	-
<i>Rhizopus</i> sp.	-	-	-
<i>Saccharomyces cerevisiae</i>	-	-	-
<i>Scopulariopsis brevicaulis</i>	-	-	-
<i>Sporothrix schenckii</i>	-	+	+
<i>Trichosporon beigelii</i>	-	-	-

+	Positive
-	Negative after 20 minute wash in 2X SSC

## 5 Table 7:

GenBank search results listing genes from other organisms having 100% identity to probes SEQ ID NO: 3 to SEQ ID NO: 23

	PROBE SEQ ID NO:	ORGANISM MATCHED	GENE MATCHED* (see note below)	ACCESSION NUMBER
<i>Aspergillus fumigatus</i>	3	-	-	-
<i>Blastomyces dermatitidis</i>	4	<i>Streptomyces verticillus</i>	bleomycin acetyl transferase	L26955
	4	<i>Giardia muris</i>	upstream of rRNA genes	X65063, S53320
	4	<i>Aspergillus nidulans</i>	uric acid-xanthine permease	X71807
	4	<i>Homo sapiens</i>	T-cell surface glycoprotein	X16996
	4	<i>Homo sapiens</i>	MIC2	M16279, M22557, J03841, M22556
<i>Candida albicans</i>	5	<i>Candida albicans</i>	28S rRNA	L28817
<i>Coccidioides immitis</i>	6	-	-	-
<i>Cryptococcus neoformans</i>	7	<i>Cryptococcus neoformans</i>	28S rRNA	L14067, L14068,
<i>Cryptococcus neoformans</i>	8	<i>Cryptococcus neoformans</i>	28S rRNA	L14067, L14068, L20964
	8	<i>Escherichia coli</i>	0111 cld	Z17241
<i>Histoplasma capsulatum</i>	9	-	-	-
<i>Aspergillus glaucus</i>	10	<i>Pseudomonas denitrificans</i>	cob genes	M62866
<i>Aspergillus niger</i>	11	-	-	-
<i>Aspergillus terreus</i>	12	Human cytomegalovirus	genome	X17403
	12	<i>Homo sapiens</i>	GABA receptor	L08485
<i>Candida glabrata</i>	13	<i>Homo sapiens</i>	Class I MHC	X03664, X03665
<i>Candida guilliermondii</i>	14	-	-	-
<i>Candida kefyr</i>	15	-	-	-

<i>Candida krusei</i>	16	<i>Pseudomonas syringae</i>	penicillin binding protein	L28837
<i>Candida lusitanae</i>	17	Chicken	AK1	D00251
	17	Mouse	IL10	M84340
<i>Candida parapsilosis</i>	18	<i>Polytomella agilis</i>	beta-2 tubulin	M33373
	18	Tobacco chloroplast	genome	Z00044, S54304
	18	<i>Aedes aegypti</i>	amylase	L03640
	18	<i>Homo sapiens</i>	chromosome 13q14	L14473
<i>Candida tropicalis</i>	19	-	-	-
<i>Pseudallescheria boydii</i>	20	<i>Drosophila melanogaster</i>	AcTr66B	X71789
		Cow	actin 2	D12816
<i>Aspergillus flavus</i>	21	-	-	-
<i>Sporothrix schenckii</i>	22	-	-	-
<i>Sporothrix schenckii</i>	23	Sulfate reducing bacteria	FMN binding protein	D21804
	23	Equine herpesvirus 1	genome	M86664

\* Note: As discussed earlier in this document, the presence of sequences similar to probes SEQ ID NO:3 to SEQ ID NO: 23 in genes not related to 28S does not have any effect on the specificity or sensitivity of our diagnostic strategy. Our species specific probes are used to analyze 28S DNA that has been previously amplified in a polymerase chain reaction with our probes SEQ ID NO: 1 and SEQ ID NO:2. These two probes will not amplify DNA from any gene other than 28S in column #4 (GENE MATCHED), and therefore no amplified DNA from these non-28S genes will be available for the hybridization of probes SEQ ID NO: 3 to SEQ ID NO: 23.

10 EXAMPLE 2. Use of method in example 1 to test clinical specimens for specific fungal organisms.

Clinical samples taken from the respiratory and gastrointestinal tract of healthy individuals almost always contain some fungal flora. Most of these fungi are non-pathogenic, but may give false positives on traditional immunochemical diagnostic tests for pathogenic fungi.

We obtained 44 clinical specimens from diverse sources ranging from sputum and incision drainage tubes, to intervertebral disc and lung biopsies. Traditional smear and culture results showed that all 44 specimens contained at least 1 type of fungus. In order to test the efficacy of

our probes, we extracted DNA from all 44 clinical samples and used probes SEQ ID NO: 1 & 2 in a polymerase chain reaction to amplify fungal 28S sequences present in these samples.

DNA was extracted from all clinical samples by our modification of the technique of Chomczynski and Sacchi which originally described the use of acid guanidinium thiocyanate-phenol-chloroform to preferentially extract RNA from cells and tissues. We replaced room temperature cell lysis by boiling lysis, and acid guanidinium thiocyanate-phenol-chloroform extraction by alkaline phenol-guanidine thiocyanate to preferentially extract DNA from cells. 1.5 ml Sarsted (Newton, North Carolina) polypropylene screw cap tubes with o-ring seals were used for the extractions. 200 ul of specimen was added to 500 ul of GPT reagent (6 M guanidine thiocyanate dissolved in 50 mM tris pH 8.3 mixed with an equal volume of phenol buffered in tris pH 8.0). This was mixed by vortexing and immediately placed in a boiling water bath for 15 minutes. The tubes were spun in a microcentrifuge for 5 seconds and 250 ul of chloroform/iso-amyl alcohol (24:1 by volume) was added and mixed by vortexing. The liquid phases were separated by centrifugation for 10 minutes and 450 ul of aqueous (upper) phase was transferred to a fresh tube. The aqueous phase was mixed with 500 ul of 100% isopropanol and placed at -20°C for at least 1 hour. At the end of this period the tubes were centrifuged for 15 minutes and the supernatant removed without disturbing the nucleic acid pellet. The pellet was washed with 500 ul of ice-cold 70% ethanol to remove traces of GPT reagent by gently inverting 2 times and then centrifuged for 5 minutes. The ethanol was removed and the pellet dried in a speed vac for 10 minutes. The pellet was resuspended in 25 ul of sterile deionized water and 5 ul was used in a 50 ul PCR amplification. The PCR was carried out as a hot-start reaction using the thermal cycling conditions for probes SEQ ID NO: 1 and SEQ ID NO: 2 described in example 1. Gel electrophoresis showed that probes SEQ ID NO: 1 and SEQ ID NO: 2 successfully amplified DNA from all 44 specimens.

The amplified DNA from each specimen was transferred to a positively charged polysulphone based membrane. We radioactively labeled our species specific probes SEQ ID

NO: 3, SEQ ID NO: 5, SEQ ID NO: 6 and SEQ ID NO: 7, and sequentially probed the membrane to test for the presence of 28S rDNA from *Aspergillus fumigatus*, *Candida albicans*, *Coccidioides immitis* and *Cryptococcus neoformans* respectively. Membrane blocking, probe hybridization and washes were done exactly as described in example 1. The results are shown in Table 8.

5

No false positives were observed, indicating a specificity of 100% for these 4 probes in the clinical specimens tested. 10 out of 12 culture positive samples for *Aspergillus fumigatus*, and 11 out of 13 samples of *Candida albicans* were identified, indicating a detection sensitivity of about 85% for these two probes. Additionally, two out of two *Coccidioides immitis* and two out of two *Cryptococcus neoformans* were correctly identified (detection sensitivity of 100%). As seen by these results, the probes described in this invention can be used on a diverse variety of clinical specimens with excellent efficacy.

10

Table 8.

15

Detection of *Aspergillus fumigatus*, *Candida albicans*, *Coccidioides immitis* and *Cryptococcus neoformans* in clinical specimens using species specific probes.

Specimen type	Smear and culture results	PCR with SEQ ID: 1, 2	SEQ ID: 3	SEQ ID: 5	SEQ ID: 6	SEQ ID: 7
U035 sputum	A. flavus	+	-	-	-	-
U069 pleura	A. fumigatus	+	+	-	-	-
U070 bronchial wash	A. flavus	+	-	-	-	-
M019 bronchial wash	A. fumigatus	+	+	-	-	-
M020 sputum	mixed fungal flora	+	-	+	-	-
X35254 sputum	C. albicans	+	-	+	-	-
M20910 sputum	A. fumigatus	+	+	-	-	-
M055 sputum	C. albicans	+	-	+	-	-
M056 abdominal	mixed fungal flora	+	-	-	-	-
M057 drainage tube	C. albicans	+	-	(-)	-	-
M059 ind. sputum	C. albicans	+	-	+	-	-
M060 ind. sputum	mixed fungal flora	+	-	-	-	-
M083 bronchial wash	C. albicans	+	-	+	-	-
M084 sputum	A. fumigatus	+	(-)	-	-	-

M085 throat	C. albicans	+	-	(-)	-	-
A001 sputum	A. fumigatus	+	(-)	-	-	-
A002 leg	Blastomyces	+	-	-	-	-
A003 leg	Blastomyces	+	-	-	-	-
A005 disc	A. fumigatus	+	+	-	-	-
A037 disc	A. fumigatus	+	+	-	-	-
A039 trachea	C. albicans	+	-	+	-	-
A040 trachea	C. albicans	+	-	+	-	-
A102 empyema	A. fumigatus	+	+	-	-	-
Y004 sputum	C. albicans	+	-	+	-	-
Y016 induced sputum	Coccidioides	+	-	-	+	-
Y028 sputum	Coccidioides	+	-	-	+	-
J003 chest	Aspergillus sp.	+	-	-	-	-
J045 bronchial wash	C. albicans	+	-	+	-	-
J046 ethmoid	yeast	+	-	-	-	-
J047 chest	A. fumigatus	+	+	-	-	-
J048 sputum	C. albicans	+	-	+	-	-
J073 lung	Aspergillus sp.	+	-	-	-	-
J074 lung	A. fumigatus	+	+	-	-	-
U017 lip	A. fumigatus	+	+	-	-	-
U033 sputum	mixed fungal flora	+	-	-	-	-
U071 sputum	C. albicans	+	-	+	-	-
U072 BA lavage	Sporothrix	+	-	-	-	-
U073 knee	Histoplasma	+	-	-	-	-
U074 mandible	Cryptococcus	+	-	-	-	+
U075 CSF	Cryptococcus	+	-	-	-	+
U076 knee	Histoplasma	+	-	-	-	-
U077 soft tissue	Histoplasma	+	-	-	-	-
U051 buccal	A. fumigatus	+	+	-	-	-
Y055 sputum	mixed fungal flora	+	-	-	-	-
+ Positive    - Negative    (-) Missed						



**EXAMPLE 3. DNA sequence based identification of unknown fungal organisms.**

Another utility of our probes is in the rapid DNA sequence based identification of a pure culture of fungus. Probes SEQ ID NO: 1 and SEQ ID NO: 2 are used in a  
5 polymerase chain reaction to amplify 28S rDNA from an unknown fungus. Probes SEQ ID NO: 1 or SEQ ID NO: 2 are then used as sequencing primers to obtain DNA sequence from this amplified 28S DNA belonging to the unknown fungus. This DNA sequence is compared to the fungal 28S DNA sequences in our database, and a sequence match at, or overlapping any one of the probe sequences in SEQ ID NO: 3 to SEQ ID NO: 74 will  
10 confirm the identity of the fungus. This technique cannot be used directly on clinical samples, as these usually contain DNA from more than one fungus, and the DNA sequence generated will consist of overlapping sequences of several organisms. This technique has utility in rapidly and reliably identifying colonies of a single fungus on culture plates, clinical specimens, food, pharmaceutical, environmental or other samples  
15 containing only one species of fungus.

**EXAMPLE 4. Capture and identification of target DNA or RNA**

All primers and probes described in this invention disclosure may be labeled with  
20 any detectable reporter or signal moiety including, but not limited to radioisotopes, enzymes, antigens, antibodies, chemiluminescent reagents and fluorescent chemicals. Additionally, these probes may be modified without changing the substance of their purpose by terminal addition of nucleotides designed to incorporate restriction sites or other useful sequences. These probes may also be modified by the addition of a capture  
25 moiety (including, but not limited to para-magnetic particles, biotin, fluorescein, dioxigenin, antigens, antibodies) or attached to the walls of microtiter trays to assist in the solid phase capture and purification of these probes and any DNA or RNA hybridized to these probes. Fluorescein may be used as a signal moiety as well as a capture moiety, the latter by interacting with an anti-fluorescein antibody.

30

A typical utility of these modifications would be as follows. Primers SEQ ID NO: 1 and SEQ ID NO: 2 would be utilized to amplify 28S rDNA from a sample, if present, as described previously. Primers would be modified so as to contain a biotin moiety at their 5' ends. A streptavidin solid phase, such as a paramagnetic particle, would be used to  
5 separate PCR products, if present, from the reaction mixture. The amplified target may be subsequently hybridized to a third probe ((SEQ ID NO: 3) to (SEQ ID NO: 74) or their complements) attached to a detectable moiety to determine which species of fungus is present in the given sample. Multiple probes, each labeled with a different detectable moiety may be used at one time to analyze the amplified target.

10

Alternatively, Primers SEQ ID NO: 1 and SEQ ID NO: 2 would be utilized to amplify 28S rDNA from a sample, if present, as above. In a separate reaction, individually, either SEQ ID NO: 1 or SEQ ID NO: 2 would be modified by attachment to a solid phase capture moiety, such as a paramagnetic particle, and SEQ ID NO: 3 to SEQ  
15 ID NO: 74 (or their complements) would be modified by addition of a detectable moiety. Alternately, in the amplicon, any sequences delimited by SEQ ID NO: 1 and SEQ ID NO: 2, including but not limited to SEQ ID NO: 3 to SEQ ID NO: 74, may be used in the design of a capture probe. One of the probes attached to a solid phase (SEQ ID NO: 1 and SEQ ID NO: 2) or any other appropriately designed sequences and one of the probes  
20 modified by attachment to a detectable moiety (SEQ ID NO: 3 to SEQ ID NO: 74 or their complements) would be hybridized together, in solution, to products of the PCR, if they had been generated. The hybrids, if present, would be captured from the solution, and analyzed by a method appropriate to the detection moiety. Detection of the hybridized probe would indicate which species of fungus was present in the given sample. Multiple  
25 probes, each labeled with a different detectable moiety may be used at one time to analyze the amplified target.

30

**EXAMPLE 5. Species-specific amplification of fungal DNA**

Another utility of the probes described in this invention is their usage as primers in the direct detection of a specific fungal species by virtue of a nucleic acid amplification  
5 reaction. In this embodiment, one primer is a universal one, such as (SEQ ID NO:1) or (SEQ ID NO:2), and the other is a species-specific primer selected from the group consisting of (SEQ ID NO:3) to (SEQ ID NO: 23) or the complements thereof. One variation of this approach is the substitution of (SEQ ID NO:1) or (SEQ ID NO:2) with any functional sequence located in proximity to the species-specific primer. Another  
10 variation of this approach is the selection of any appropriate species specific primer pair from SEQ ID NO: 24 to SEQ ID NO: 74.

## SEQUENCE LISTING

## (1) GENERAL INFORMATION:

## 5 (i) APPLICANT:

(A) Sandhu, Gurpreet S.

(B) Kline, Bruce C.

## (ii) TITLE OF INVENTION:

10 Nucleic Acid Probes for the Detection and Identification of Fungi

## (iii) NUMBER OF SEQUENCES: 23

## (iv) CORRESPONDENCE ADDRESS:

15 (A) ADDRESSEE: Ciba Corning Diagnostics Corp.

(B) STREET: 63 North Street

(C) CITY: Medfield

(D) STATE: Massachusetts

(E) COUNTRY: USA

20 (F) ZIP: 02052

## (v) COMPUTER READABLE FORM:

(A) MEDIUM TYPE: Diskette 3.5 inch, 1.44 Mb storage

(B) COMPUTER: IBM PS/2

25 (C) OPERATING SYSTEM: MS-DOS 6.2

(D) SOFTWARE: Word 6.0

## (vi) CURRENT APPLICATION DATA:

(A) APPLICATION NUMBER:

30 (B) FILING DATE:

## (C) CLASSIFICATION:

## (vii) PRIOR APPLICATION DATA:

## 5 (viii) ATTORNEY INFORMATION:

(A) NAME: Morgenstern, Arthur S.

(B) REGISTRATION NUMBER: 28,244

(C) DOCKET NUMBER: CCD-180

## 10 (ix) TELECOMMUNICATION INFORMATION:

(A) TELEPHONE: 508 359-3836

(B) TELEFAX: 508 359-3885

## (2) INFORMATION FOR SEQ ID NO 1:

15

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 20

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

20

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for fungal organisms

(iii) HYPOTHETICAL: No

25

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 1:

30 GTGAAATTGT TGAAAGGGAA

20

## (3) INFORMATION FOR SEQ ID NO 2:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 18  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: Nucleic acid probe for fungal organisms

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 2:

GACTCCTTGG TCCGTGTT 18

## (4) INFORMATION FOR SEQ ID NO 3:

20

## (i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 14  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
25 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Aspergillus fumigatus*

(iii) HYPOTHETICAL: No

30

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 3:

5 CTCGGAATGT ATCA 14

(5) INFORMATION FOR SEQ ID NO 4:

(i) SEQUENCE CHARACTERISTICS:

10 (A) LENGTH: 13  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

15 (ii) MOLECULE TYPE: Nucleic acid probe for *Blastomyces dermatitidis*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

20

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 4:

ACTCCCCAC GGG 13

25 (6) INFORMATION FOR SEQ ID NO 5:

(i) SEQUENCE CHARACTERISTICS:

30 (A) LENGTH: 14  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida albicans*

5 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 5:

10

CCTCTGACGA TGCT

14

(7) INFORMATION FOR SEQ ID NO 6:

15 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

20

(ii) MOLECULE TYPE: Nucleic acid probe for *Coccidioides immitis*

(iii) HYPOTHETICAL: No

25 (iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 6:

TCTGGCGGTT GGTT

14

30 (8) INFORMATION FOR SEQ ID NO 7:



## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

5 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Cryptococcus neoformans*

(iii) HYPOTHETICAL: No

10

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 7:

15 CTCCTGTCGC ATAC

14

## (9) INFORMATION FOR SEQ ID NO 8:

## (i) SEQUENCE CHARACTERISTICS:

20 (A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

25 (ii) MOLECULE TYPE: Nucleic acid probe for *Cryptococcus neoformans*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

30 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 8:

AGTTCTGATC GGTG

14

(10) INFORMATION FOR SEQ ID NO 9:

5

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

10

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Histoplasma capsulatum*

(iii) HYPOTHETICAL: No

15

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 9:

20 CAATCCCCCG CGGC

14

(11) INFORMATION FOR SEQ ID NO 10:

(i) SEQUENCE CHARACTERISTICS:

25

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

30

(ii) MOLECULE TYPE: Nucleic acid probe for *Aspergillus glaucus*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

5

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 10:

GTGTCATGCG GCCA

14

10 (12) INFORMATION FOR SEQ ID NO 11:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

15

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Aspergillus niger*

20

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 11:

25

CCCTGGAATG TAGT

14

(13) INFORMATION FOR SEQ ID NO 12:

(i) SEQUENCE CHARACTERISTICS:

30

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

5 (ii) MOLECULE TYPE: Nucleic acid probe for *Aspergillus terreus*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

10

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 12:

GCTTCGGCCC GGTG

14

15 (14) INFORMATION FOR SEQ ID NO 13:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

20

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida glabrata*

25

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 13:

30 CTTGGGACTC TCGC

14

## (15) INFORMATION FOR SEQ ID NO 14:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 14  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: Nucleic acid probe for *Candida guilliermondii*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 14:

ATATTTTGTG AGCC

14

20 (16) INFORMATION FOR SEQ ID NO 15:

## (i) SEQUENCE CHARACTERISTICS:

- 25 (A) LENGTH: 14  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida kefyr*

30 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 15:

5

TTCGGCTTTC GCTG

14

(17) INFORMATION FOR SEQ ID NO 16:

10

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

15

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida krusei*

(iii) HYPOTHETICAL: No

20

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 16:

GGGATTGCGC ACCG

14

25 (18) INFORMATION FOR SEQ ID NO 17:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

30

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida lusitaniae*

5 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 17:

10

GCCTCCATCC CTTT

14

(19) INFORMATION FOR SEQ ID NO 18:

15 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

20

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida parapsilosis*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 18:

ATAAGTGCAA AGAA

14

30 (20) INFORMATION FOR SEQ ID NO 19:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

5 (C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Candida tropicalis*

10 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 19:

15

AGAATTGCGT TGGA

14

## (21) INFORMATION FOR SEQ ID NO 20:

20 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

25 (ii) MOLECULE TYPE: Nucleic acid probe for *Pseudallescheria boydii*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

30



(v) SEQUENCE DESCRIPTION: SEQ ID NO: 20:

GCGATGGGAA TGTG

14

5 (22) INFORMATION FOR SEQ ID NO 21:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

10 (C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Aspergillus flavus*

15 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 21:

20

AGACTCGCCT CCAG

14

(23) INFORMATION FOR SEQ ID NO 22:

(i) SEQUENCE CHARACTERISTICS:

25 (A) LENGTH: 14

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

30 (ii) MOLECULE TYPE: Nucleic acid probe for *Sporothrix schenckii*

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

5

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 22:

CGGACCACCC GGCG

14

10 (24) INFORMATION FOR SEQ ID NO 23:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 14

(B) TYPE: nucleic acid

15

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: Nucleic acid probe for *Sporothrix schenckii*

20

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 23:

CGGCGGCATG CCCC

14

25

(25) INFORMATION FOR SEQ ID NO 24:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 208

30

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Acremonium* species specific region of 28S gene.

5

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

10

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 24:

15

GACCAGACTT GGGCTCGGTG AATCATCCGG CGTTCTCGCC GGTGCACTTT  
GCCGTCCCAG GCCAGCATCA GTTCGCGCCG GGGGATAAAG GTTTCGGGAA  
TG TAGCTCCT TCGGGAGTGT TATAGCCCGT TCGGTAATAC CCTGGCGTGG  
ACTGAGGTCC GCGCTCTGCA AGGATGCTGG CGTAATGGTC ATCAGTGACC  
CGTCTTGA

20 (26) INFORMATION FOR SEQ ID NO 25:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

(B) TYPE: nucleic acid

25

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus clavatus* specific region of 28S gene.

30

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 25:

5  
 GACCAGACTC GCTCGCGGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
 CCGTGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCTCCGGAA  
 TGTATCACCT CTCGGGGTGT CTTATAGCCG GGGGTGCAAT GCGGCCTGCC  
 TGGACCGAGG AACGCGCTTC GGCTCGGACG CTGGCGTAAT GGTCGTAAAT  
 10 GACCCGTCTT GA

(27) INFORMATION FOR SEQ ID NO 26:

15 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

20 (ii) MOLECULE TYPE: *Aspergillus flavus* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 26:

30 GACCAGACTC GCCTCCAGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
 CTGGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCTCCCGGAA  
 TG TAGTGCCC TYCGGGGCAC CTTATAGCCG GGAGTGCAAT GCGGCCAGCC  
 TGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCATAAT GGTCGYAAAC  
 GACCCGTCTT GA

## (28) INFORMATION FOR SEQ ID NO 27:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

5

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus fumigatus* specific region of 28S gene.

10

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 27:

20 GACCAGACTC GCCCGCGGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
CCGTGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCTCGGAA  
TGTATCACCT CTCGGGGTGT CTTATAGCCG AGGGTGCAAT GCGGCCTGCC  
TGGACCGAGG AACGCGCTTC GGCTCGGACG CTGGCGTAAT GGTCGTAAAT  
GACCCGTCTT GA

## 25 (29) INFORMATION FOR SEQ ID NO 28:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

30

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus glaucus* specific region of 28S gene.

(iii) HYPOTHETICAL: No

5 (iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 28:

10 GACCAGACTC GCTTCGGGGG TTCAGCCGGC TTTCGGGCGG GTGTACTTCC  
CCGGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCCTGGAA  
TGTAACGCCT CTCGGGGCGC CTTATAGCCA GGGGTGTCAT GCGGCCAGCC  
TGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCATAAT GGTCTGTAAC  
GACCCGTCTT GA

15

(30) INFORMATION FOR SEQ ID NO 29:

(i) SEQUENCE CHARACTERISTICS:

- 20 (A) LENGTH: 213  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

25 (ii) MOLECULE TYPE: *Aspergillus nidulans* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

30

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 29:

5 GACCAGACTC GGCCCCGGGG TTCARCCAGC ACTCGTGCTG GTGTACTTCC  
CCGGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCCAGGAA  
TGTATCGCCC TCCGGGGTTG TCTTATAGCC TGGGGTGCAA TGCGGCCAGC  
CCGGACCGAG GAACGCGCTT CGGCACGGAC GCTGGCGTAA TGGTCGCAAA  
CGACCCGTCT TGA

10 (31) INFORMATION FOR SEQ ID NO 30:

(i) SEQUENCE CHARACTERISTICS:

- 15 (A) LENGTH: 212  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus niger* specific region of 28S gene.

20 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 30:

25 GACCAGACTC GCCC GCGGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
CCGTGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCCAGGAA  
TG TAGTRCCC TCCGGGGYAC CTTATAGCCA GGGGTGCAAT GCGGCCAGCC  
TGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCATAAT GGTGCTAAAC  
30 GACCCGTCTT GA

## (32) INFORMATION FOR SEQ ID NO 31:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

5 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus ochraceus* specific region of 28S gene.

10

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 31:

15

GACCAGACTC GCCCGCGGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
CCGCGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCGCCCGGAA  
TGTCAGCACC TTCGGGGTGC CTTATAGCCG GGGGTGCAAT GCGGCCAGCC  
20 TGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCATAAT GGTCGTAAAC  
GACCCGTCTT GA

## (33) INFORMATION FOR SEQ ID NO 32:

25

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH:

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

30 (D) TOPOLOGY: linear



(ii) MOLECULE TYPE: *Aspergillus terreus* specific region of 28S gene.

(iii) HYPOTHETICAL: No

5 (iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 32:

10 AACCAGACTC GCTCGCGGGG TTCAGCCGGG CTCGGCCCCG GTGTACTTCC  
CCGCGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCTCCGGAA  
TGTAGCGCCC TTCGGGGCGC CTTATAGCCG GGGGTGCAAT GCGGCCAGCC  
TGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCATAAT GGTTGTAAAC  
GACCCGTCTT GA

15

(34) INFORMATION FOR SEQ ID NO 33:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 213

20 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus unguis* specific region of 28S gene.

25

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

30 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 33:

GACCAGACTC GGCCTCGGGG TTCAGCCAGC ACTCGTGCTG GTGTACTTCC  
 CCGGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCCAGGAA  
 TGTATCACCC TCCGGGGTTG TCTTATAGCC TGGGGTGCAA TGCGGCCAGC  
 5 CTGGACCGAG GAACGCGCTT CGGCACGGAC GCTGGCATAA TGGTTGCAAA  
 CGACCCGTCT TGA

(35) INFORMATION FOR SEQ ID NO 34:

10

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 212

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

15

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Aspergillus ustus* specific region of 28S gene.

(iii) HYPOTHETICAL: No

20

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 34:

25

GACCAGACTC GGCCCCGGGG TTCAGCCAGC ACTCGTGCTG GTGTACTTCC  
 CCGGGGGCGG GCCAGCGTCG GTTTGGGCGG CCGGTCAAAG GCCCCAGGAA  
 TGTGTCGCCC TCCGGGGCGT CTTATAGCCT GGGGTGCAAT GCGGCCAGCC  
 CGGACCGAGG AACGCGCTTC GGCACGGACG CTGGCGTAAT GGTGCGAAAC  
 GACCCGTCTT GA

30

## (36) INFORMATION FOR SEQ ID NO 35:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 208

5 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Beauveria* species specific region of 28S gene.

10

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 35:

20 GACCAGACTT GGGCTTGGTT GATCATCCGG GGTTCCTCCCC GGTGCACTCT  
TCCGGCCCAG GCCAGCATCA GTTCGCCCTG GGGGACAAAG GCTTCGGGAA  
CGTGGCTCTC TCCGGGGAGT GTTATAGCCC GTTGCGTAAT ACCCTGTGGC  
GGACTGAGGT TCGCGCATTG GCAAGGATGC TGGCGTAATG GTCATCAGTG  
ACCCGTCT

## 25 (37) INFORMATION FOR SEQ ID NO 36:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 213

30 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Bipolaris* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

5

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 36:

10

AGCCAGACTT GCTTGCAGTT GCTCATCCGG GCTTTTGCCC  
GGTGCACCTCT TCTGCAGGCA GGCCAGCATC AGTTTGGGCG  
GTGGGATAAA GGTCTCTGTC ACGTACCTTC CTTCGGGTTG  
GCCATATAGG GGAGACGTCA TACCACCAGC CTGGACTGAG  
15 GTCCGCGCAT CTGCTAGGAT GCTGGCGTAA TGGCTGTAAG  
CGGCCCGTCT TGA

(38) INFORMATION FOR SEQ ID NO 37:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 105

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

25

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Blastoschizomyces* species specific region of 28S gene.

30

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 37:

5       TGAAATTGTT GAAAGGGAAG GCGATGGTAG GAATAAGAGG CTGCGGTTTG  
AAATAATTGT TTTTCGGGCC ACGGTCTCCT GAGCCTGCTT TCGCACCCGT  
CTTGA

## 10   (39) INFORMATION FOR SEQ ID NO 38:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 214

(B) TYPE: nucleic acid

15       (C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Blastomyces dermatitidis* specific region of 28S  
gene.

20

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

## 25       (v) SEQUENCE DESCRIPTION: SEQ ID NO: 38:

30       GACCAGAGTC GGCCGTGGGG GTTCAGCGGG CATTCGTTGC CCGTGCACTC  
CCCCACGGGC GGGCCAGCGT CGGTTTCGAC GGCCGGTCAA AGGCCCCCGG  
AATGTGTCGC CTCTCGGGGC GTCTTATAGC CGGGGGTGCA ATGCGGCCAG  
TCGGGACCGA GGAACGCGCT TCGGCACGGA CGCTGGCTTA ATGGTCGTAA  
GCGACCCGTC TTGA

## (40) INFORMATION FOR SEQ ID NO 39:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 213  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: *Chrysosporium* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 39:

20 AACCAGACTT GCGCGCGGCC GATCATCCGG TGTTCTCACC GGTGCACTCG  
GCCGTGCTCA GGCCAGCATC GGTTTGGCG GCTGGATAAA GGCCCTAGGA  
ATGTGGCTCC TCTCGGGGAG TGTTATAGCC TAGGGTGCAA TGCAGCCTGC  
TGGGACCGAG GACCGCGCTT CGGCTAGGAT GCTGGCGTAA TGGTTGTAAG  
CGGCCCGTCT TGA

25

## (41) INFORMATION FOR SEQ ID NO 40:

## (i) SEQUENCE CHARACTERISTICS:

- 30 (A) LENGTH: 207  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Cladosporium* species specific region of 28S gene.

5 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 40:

10

AACCAGACTT GCTCGCGGTG TTCCGCCGGT CTTCTGACCG GTCTACTCGC  
CGCGTTGCAG GCCAGCATCG TCTGGTGCCG CTGGATAAGA CTTGAGGAAT  
GTAGCTCCCT CGGGAGTGTT ATAGCCTCTT GTGATGCAGC GAGCGCCGGG  
15 CGAGGTCCGC GCTTCGGCTA GGATGCTGGC GTAATGGTCG TAATCCGCCC  
GTCTTGA

(42) INFORMATION FOR SEQ ID NO 41:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 213

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

25

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Curvularia* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

30

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 41:

5 AGCCAGACTT GCTTGCACTT GCTCATCCGG GCTTTTGCCC GGTGCACTCT  
TCTGCAGGCA GGCCAGCATC AGTTTGGGCG GTGGGATAAA GGTCTCTGAC  
ACGTTCCCTTC CTTCCGGGTTG GCCATATAGG GGAGACGTCA TACCACCAGC  
CTGGACTGAG GTCCGCGCAT CTGCTAGGAT GCTGGCGTAA TGGCTGTAAG  
10 CGGCCCCGTCT TGA

## (43) INFORMATION FOR SEQ ID NO 42:

## (i) SEQUENCE CHARACTERISTICS:

15 (A) LENGTH: 213  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

20 (ii) MOLECULE TYPE: *Candida albicans* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 42:

30 GATCAGACTT GGTATTTTGC ATGCTGCTCT CTCGGGGGCG GCCGCTGCGG  
TTTACCGGGC CAGCATCGGT TTGGAGCGGC AGGATAATGG CGGAGGAATG  
TGGCACGGCT TCTGCTGTGT GTTATAGCCT CTGACGATGC TGCCAGCCTA  
GACCGAGGAC TGCGGTTTTT AACCTAGGAT GTTGGCATAA TGATCTTAAG



TCGCCCCGTCT TGA

## (44) INFORMATION FOR SEQ ID NO 43:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 223  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

- 10 (ii) MOLECULE TYPE: *Candida glabrata* specific region of 28S gene.

## (iii) HYPOTHETICAL: No

## (iv) ANTISENSE: No

15

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 43:

20 GATCAGACAT GGTGTTTTGC GCCCCTTGCC TCTCGTGGGC TTGGGACTCT  
CGCAGCTCAC TGGGCCAGCA TCGGTTTTGG CGGCCGAAA AAACCTAGGG  
AATGTGGCTC TGCGCCTCGG TGTAGAGTGT TATAGCCCTG GGGAATACGG  
CCAGCCGGGA CCGAGGACTG CGATACTTGT TATCTAGGAT GCTGGCATAA  
TGGTTATATG CCGCCCCGTCT TGA

25

## (45) INFORMATION FOR SEQ ID NO 44:

## (i) SEQUENCE CHARACTERISTICS:

- 30 (A) LENGTH: 212  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Candida guilliermondii* specific region of 28S gene.

5 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 44:

10

GATCAGACTC GATATTTTGT GAGCCTTGCC TTCGTGGCGG GGTGACCCGC  
AGCTTATCGG GCCAGCATCG GTTTGGGCGG TAGGATAATG GCGTAGGAAT  
GTGACTTTTC TCGGTGAAG TGTATAGCC TGCCTTGATG CTGCCTGCCT  
15 AGACCGAGGA CTGCGATTTT ATCAAGGATG CTGGCATAAT GATCCCAAAC  
CGCCCGTCTT GA

(46) INFORMATION FOR SEQ ID NO 45:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 214

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

25

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Coccidioides immitis* specific region of 28S gene.

(iii) HYPOTHETICAL: No

30

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 45:

5       AACCAGACTC GGTCTGTTGGG GCTCAGCGGG CATGAGTGCC CGTGTACTCC  
       CCCATGCTCC GGGCCAGCAT CAGTTCTGGC GGTGTTTAA AGGCCTCTGG  
       AATGTATCGT CCTCCGGGAC GTCTTATAGC CAGGGGCGCA ATGCGGCCAG  
       CCGGGACTGA GGAACGCGCT TCGGCACGGA TGCTGGCATA ATGGTTGTAA  
       GCGGCCCGTC TTGA

10

## (47) INFORMATION FOR SEQ ID NO 46:

## (i) SEQUENCE CHARACTERISTICS:

- 15       (A) LENGTH: 187  
       (B) TYPE: nucleic acid  
       (C) STRANDEDNESS: single  
       (D) TOPOLOGY: linear

20       (ii) MOLECULE TYPE: *Candida kefyr* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 46:

30       GATCAGACAT GGCCTTTGCT TCGGCTTTCTG CTGGGCCAGC ATCAGTTTAA  
       GCGGTTGGAT AAATCCTCGG GAATGTGGCT CTGCTTCGGT AGAGTGTAT  
       AGCCCGTGGG AATACAGCCA GCTGGGACTG AGGATTGCGA CTTTGTCAA  
       GGATGCTGGC GTAATGGTTA AATGCCGCCG GTCTTGA

## (48) INFORMATION FOR SEQ ID NO 47:

## 5 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 213

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

10

(ii) MOLECULE TYPE: *Candida krusei* specific region of 28S gene.

(iii) HYPOTHETICAL: No

15 (iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 47:

20 CGCCCGACAT GGGGATTGCG CACCGCTGCC TCTCGTGGGC GGCCTCTGG  
GCTTTCCTG GGCAGCATC GGTTCCTGCT GCAGGAGAAG GGGTCTGGA  
ACGTGGCTCT TCGGAGTGTT ATAGCCAGGG CCAGATGCTG CGTGCGGGGA  
CCGAGGACTG CGGCCGTGTA GGTCACGGAT GCTGGCAGAA CGGCGCAACA  
25 CCGCCCGTCT TGA

25

## (49) INFORMATION FOR SEQ ID NO 48:

## (i) SEQUENCE CHARACTERISTICS:

30 (A) LENGTH: 236

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Cryptococcus laurentii* specific region of 28S gene.

5

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

10

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 48:

15

AGTCAGTCGT GTCTGGGAGG CTCAGCCGGT TCTGCCGGTG TATTCCTCTC  
AGACGGGTCA ACATCAGTTT TGTCCGACGG ATAATGGCGG CGGGAAAGTA  
GCACCTCCGG GTGTGTTATA GCCCGCTGTC GCATACGCCG GATGAGACTG  
AGGCATGCAG CTCGCCTTTA TGGCAGGGGT TCGCCCACTT TCGAGCTTAG  
GATGTTGACG TAATGGCTTT AAACGACCCG TCTTGA

20 (50) INFORMATION FOR SEQ ID NO 49:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 173

(B) TYPE: nucleic acid

25

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Candida lusitanae* specific region of 28S gene.

30

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 49:

5

AAGCAGACAC GGTTTTACCG GGCCAGCGTC GAAAAGGGGG GAGGAACAAG  
AACTCGAGAA TGTGGCGCGC ACCTTCGGGY GCGCGTGTTA TAGCTCGTGT  
TGACGCCTCC ATCCCTTTTC GAGGCCTGCG ATTCTAGGAC GCTGGCGTAA  
TGGTTGCAAG CCGCCCGTCT TGA

10

(51) INFORMATION FOR SEQ ID NO 50:

(i) SEQUENCE CHARACTERISTICS:

15

- (A) LENGTH: 238
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

20

(ii) MOLECULE TYPE: *Cryptococcus neoformans* var *gattii* (serotype B)  
specific region of 28S gene.

(iii) HYPOTHETICAL: No

25

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 50:

30

AGTCAGTCGT GTCTATTGGG TTCAGCCAGC TCTGCTGGTG TATTCCCTTT  
AGACGGGTCA ACATCAGTTC TGATCGGTGG ATAAGGGCTG GAGGAATGTG  
GCACTCTTCG GGGTGTGTTA TAGCCTCCTG TCGCATAAC TGGTTGGGAC  
TGAGGAATGC AGCTCGCCTT TATGGCCGGG GTTCGCCCAC GTTCGAGCTT

AGGATGTTGA CAAAATGGCT TTAAACGACC CGTCTTGA

(52) INFORMATION FOR SEQ ID NO 51:

5 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 238

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

10

(ii) MOLECULE TYPE: *Cryptococcus neoformans* (serotype A) specific region of 28S gene.

(iii) HYPOTHETICAL: No

15

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 51:

20

AGTCAGTCGT GTCTATTGGG TTCAGCCAGT TCTGCTGGTG TATTCCCTTT  
AGACGGGTCA ACATCAGTTC TGATCGGTGG ATAAGGGCTG GGGGAATGTA  
GCACTCTTCG GAGTGTGTTA TAGCCTCCTG TCGCATAAC TGGTTGGGAC  
TGAGGAATGC AGCTCGCCTT TATGGCCGGG GTTCGCCCAC GTTCGAGCTT  
25 AGGATGTTGA CAAAATGGCT TTAAACGACC CGTCTTGA

(53) INFORMATION FOR SEQ ID NO 52:

30 (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 211

- (B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

5 (ii) MOLECULE TYPE: *Candida parapsilosis* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

10

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 52:

15 GATCAGACTT GGTATTTTGT ATGTTACTCT CTCGGGGGTG GCCTCTACAG  
TTTACCGGGC CAGCATCAGT TTGAGCGGTA GGATAAGTGC AAAGAAATGT  
GGCACTGCTT CGGTAGTGTG TTATAGTCTT TGTCGATACT GCCAGCTTAG  
ACTGAGGACT GCGGCTTCGG CCTAGGATGT TGGCATAATG ATCTTAAGTC  
GCCCCGTCTTG A

20 (54) INFORMATION FOR SEQ ID NO 53:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 238  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

25

(ii) MOLECULE TYPE: *Cryptococcus terreus* specific region of 28S gene.

30

(iii) HYPOTHETICAL: No



(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 53:

5 AGTCAGTCAT GTCTATTGGA CTCAGCCGGT TCTGCCGGTG TACTTCCTTT  
AGATGGGGTC AACATCAGTT TTGATCGCTG GAAAAGGGCA GGAGGAATGT  
AGCACTCTCG GGTGAACTTA TAGCCTTCTG TCGTATACAG TGGTTGGGAC  
TGAGGAACGC AGCATGCCTT TATGGCCGGG GTTCGCCCAC GTACATGCTT  
10 AGGATGTTGA CATAATGGCT TTAAACGACC CGTCTTGA

(55) INFORMATION FOR SEQ ID NO 54:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 211

15 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Candida tropicalis* specific region of 28S gene.

20

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 54:

GATCAGACTT GGTATTTTGT ATGTTACTTC TTCGGGGGTG GCCTCTACAG  
TTTATCGGGC CAGCATCAGT TTGGGCGGTA GGAGAATTGC GTTGGAATGT  
GGCACGGCTT CGGTTGTGTG TTATAGCCTT CGTCGATACT GCCAGCCTAG  
30 ACTGAGGACT GCGGTTTATA CCTAGGATGT TGGCATAATG ATCTTAAGTC  
GCCCCGTCTT A

## (56) INFORMATION FOR SEQ ID NO 55:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 211

5 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: *Fusarium* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 55:

20 GACCAGACTT GGGCTTGGTT AATCATCTGG GGTTCTCYCC AGTGCACTTT  
TCCAGTCCAG GCCAGCATCA GTTTTCSCCG GGGGATAAAG RCTTCGGGAA  
TGTGGCTCYC YYCGGGGAGT GTTATAGCCC GTTGYGTAAT ACCCTGGBGG  
GGACTGAGGT TCGCGCWTCT GCAAGGATGC TGGCGTAATG GTCATCAACG  
ACCCGTCTTG A

## 25 (57) INFORMATION FOR SEQ ID NO 56:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 238

30 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Filobasidium capsuligenum* specific region of 28S gene.

5 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 56:

10

AGTCAGTCAT GTCTATTGGA CTCAGCCGGT TCTGCCGGTG TATTTCTTTT  
AGATGGGGTC AACATCAGTT TTGACCGTTG GATAAAGGCA GGAAGAATGT  
AGCACTCTCG GGTGAACTTA TAGCTTCTTG TCACATACAA TGGTTGGGAC  
15 TGAGGAACGC AGCATGCCTT TATGGCCGGG ATTCGTCCAC GTACATGCTT  
AGGATGTTGA CATAATGGCT TTAAACGACC CGTCTTGA

(58) INFORMATION FOR SEQ ID NO 57:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 238

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

25

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Filobasidiella neoformans* var *bacillispora* (serotype C) specific region of 28S gene.

30 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 57:

5

AGTCAGTCGT GTCTATTGGG TTCAGCCAGC TCTGCTGGTG TATTCCTTT  
AGACGGGTCA ACATCAGTTC TGATCGGTGG ATAAGGGCTG GAGGAATGTG  
GCACTCTTCG GGGTGTGTTA TAGCCTCCTG TCGCATAAC TGGTTGGGAC  
TGAGGAATGC AGCTCGCCTT TATGGCCGGG GTTCGCCCAC GTTCGAGCTT  
10 AGGATGTTGA CAAAATGGCT TTAAACGACC CGTCTTGA

(59) INFORMATION FOR SEQ ID NO 58:

15

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 238

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

20

(ii) MOLECULE TYPE: *Filobasidiella neoformans* var *neoformans* (serotype  
D) specific region of 28S gene.

(iii) HYPOTHETICAL: No

25

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 58:

30

AGTCAGTCGT GTCTATTGGG TTCAGCCAGT TCTGCTGGTG TATTCCTTT

AGACGGGTCA ACATCAGTTC TGATCGGTGG ATAAGGGCTG GAGGAATGTG  
 GCACTCTTCG GGGTGTGTTA TAGCCTCCTG TCGCATACAC TGGTTGGGAC  
 TGAGGAATGC AGCTCGCCTT TATGGCCGGG GTTCGCCCAC GTTCGAGCTT  
 AGGATGTTGA CAAAATGGCT TTAAACGACC CGTCTTGA

5

## (60) INFORMATION FOR SEQ ID NO 59:

## (i) SEQUENCE CHARACTERISTICS:

10

(A) LENGTH: 236

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

15

(ii) MOLECULE TYPE: *Filobasidium uniguttulatum* specific region of 28S gene.

(iii) HYPOTHETICAL: No

20

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 59:

25

AGTCAGTCGT GCTCAATGGA CTCAGCCGTT CTGCGGTGTA TTCCATTGG  
 GTGGGGTCAA CATCAGTTTT GATCGCTGGA TAAAGGCAGG AGGAATGTAG  
 CACCCCCGGG TGAACCTATA GCCTCTTGTC ACATACAGTG GTTGGGACTG  
 AGGAACGCAG CATGCCTTTA TGGCCGGGAT TCGTCCACGT ACATGCTTAG  
 GATGTTGACA TAATGGCTTT AAACGACCCG TCTTGA

30

## (61) INFORMATION FOR SEQ ID NO 60:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 204  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Geotrichum* species specific region of 28S gene.

10 (iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 60:

20 AATCAGACTT GGTGCTGTTG TTCAACTRTG TTTCGGCATA GTGTACTCAG  
CAGTACTAGG CCAAGGTGGG GTGTTTGGGA GTGAAAAAGA AGTAGGAACG  
TAACTCTTCG GAGTGTTATA GCCTACTTTC ATAGCTCCTC AGGCGCCTCA  
GGACTGCGCT TCGGCAAGGA CCTTGGCATA ATGATTCTAT ACCGCCCCGTC  
TTGA

25 (62) INFORMATION FOR SEQ ID NO 61:

## (i) SEQUENCE CHARACTERISTICS:

- 30 (A) LENGTH: 214  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Histoplasma capsulatum* specific region of 28S gene.

(iii) HYPOTHETICAL: No

5

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 61:

10

GAYCAGAGTC GGCCGYGGGG GTTCAGCGGG CATTGTTGC CCGTGCAATC  
CCCCGCGGCC GGGCCAGCGT CGGTTTCGAC GGCCGGTCAA AGGCCCCCGG  
AATGTGTGCG CTCTCGGGGC GTCTTATAGC CGGGGGTGCA ATGCGGCCAG  
TCGGGACCGA GGAACGCGCT CCGGCACGGA CGCTGGCTTA ATGGTCGTCA  
15 GCGACCCGTC TTGA

(63) INFORMATION FOR SEQ ID NO 62:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 215

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

25

(ii) MOLECULE TYPE: *Malbranchea* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

30

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 62:

5 AGACAGACTC GAGCGCGGGG GCTCAGCGGG TATTGTTATG CCCGTGCACT  
CCCCCGCGCC CGGGCCAGCA TCAGTTTGG CGGCCGGTCA AAGGCCCTTG  
GAATGTATCG TCCTCCGGGA CGTCTTATAG CCAAGGGTGC AATGCGGCCA  
GCCGGGACTG AGGAACGCGC TTCGGCACGG ATGCTGGCGT AATGGCTGTA  
AGCGGCCCGT CTTGA

10

## (64) INFORMATION FOR SEQ ID NO 63:

## (i) SEQUENCE CHARACTERISTICS:

15 (A) LENGTH: 237  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

20 (ii) MOLECULE TYPE: *Mucor* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

25

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 63:

30 AGCCAGACTG GTTTGACTGT AATCAACCTA GAATTCGTTT TGGGTGCACT  
TGCAGTCTAT ACCTGCCAAC AACAGTTTGA TTTGGAGGAA AAAATTAGTA  
GGAATGTAGC CTCTCGAGGT GTTATAGCCT ACTATCATAC TCTGGATTGG  
ACTGAGGAAC GCAGCGAATG CCWTTAGGCR AGATTGCTGG GTGCTTTCGC  
TAATAAATGT TAGAATTTCT GCTTCGGGTG GTGCTAA



## (65) INFORMATION FOR SEQ ID NO 64:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 209

5 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Paecilomyces* species specific region of 28S gene.

10

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 64:

20 GACCAGACTT GGGCCCGGTG GATCATCCAG CGTTCCTCGCT GGTGCACTCC  
GCCGGGTTCA GGCCAGCATC AGTTCGCCGC GGGGGAAAAA GGCTTCGGGA  
ACGTGGCTCC TACGGGAGTG TTATAGCCCG TTGCATAATA CCCTGGGGCG  
GACTGAGGTT CGCGCTCCGC AAGGATGCTG GCGTAATGGT CATCAGCGAC  
CCGTCTTGA

## 25 (66) INFORMATION FOR SEQ ID NO 65:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 199

(B) TYPE: nucleic acid

30 (C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Penicillium* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

5

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 65:

10

GACCAGACTC GCCCACGGGG TTCAGCCGGC ATTCGTGCCG GTGTACTTCC  
CCGCGGGCGG GCCAGCGTCG GTTTGGKCGG CCGGTCAAAG GCCCTCGGAA  
TRTAACGCCC CCCGGGGCGT CTTATAGCCG AGGGTGCCAT GCGGCCAGCM  
CAGACCGAGG AACGCGCTTC GGCTCGGACG CTGGCATAAT GGTTCGTAAA

15

(67) INFORMATION FOR SEQ ID NO 66:

(i) SEQUENCE CHARACTERISTICS:

20

(A) LENGTH: 210

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

25

(ii) MOLECULE TYPE: *Pseudallescheria boydii* region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

30

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 66:

5  
GACCAGACTT GTGCCGTCG AATCAGCCGC CGCTCGTCGG CGGCGCACTT  
CGGCGGGCTC AGGCCAGCAT CAGTTCGCTG CAGGGGGAGA AAGGCGATGG  
GAATGTGGCT CTTCGGAGTG TTATAGCCCG CCGCGCAATA CCCCTCGGCG  
GACTGAGGAC CGCGCATCTG CAAGGATGCT GGCCTAATGG TCGTCAGCGA  
CCCGTCTTGA

10 (68) INFORMATION FOR SEQ ID NO 67:

**(i) SEQUENCE CHARACTERISTICS:**

(A) LENGTH: 244

**(B) TYPE: nucleic acid**

15 (C) STRANDEDNESS: single

**(D) TOPOLOGY: linear**

(ii) MOLECULE TYPE: *Rhizopus* species (NO: 1) specific region of 28S gene.

20

**(iii) HYPOTHETICAL: No**

(iv) ANTISENSE: No

25. (v) SEQUENCE DESCRIPTION: SEQ ID NO: 67:

30 AGCCGACTG GCTTGTCTGT AATCAATCTA GGTTCGTGC CTGGATGCAC  
TTGCAGACTA TTTGCCTGCC AACGACAATT TTTTTGAGT GTAAAAACTA  
TTGGAAATGT GGCCAATATT TATTTATTGG TGTATAGTC CTTTAGAAAA  
TACCTTGAAT TGGATTGAGG AACGCAGCGA ATGCTTCTCT TTAGGCGAA  
AGTCCTTTTAT TGGGATTTAC GGATCAGACT GTGGCATTGT CACA

## (69) INFORMATION FOR SEQ ID NO 68:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 215

5 (B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: *Rhizopus* species (NO: 2) specific region of 28S gene.

(iii) HYPOTHETICAL: No

15 (iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 68:

20 AGCCAGACTG GCTTGTCTGT AATCAATCTA GGCTTCGGCC TGGATGCACT  
TGCAGGCTAT GCCTGCCAAC GACAATTTGA CTTGAGGGAA AAAACTAGGG  
GAAATGTGGC CCACTTGTGG GTGTTATAGT CCCTTAGAAA ATACCTTGGG  
TTGGATTGAG GAACGCAGCG AATGCTTATT GGCGAGTTTT CCAGGAAGGT  
TTTCTGAGGT ACTAC

## (70) INFORMATION FOR SEQ ID NO 69:

25

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 215

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

30 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Rhizopus* species (NO: 3) specific region of 28S gene.

(iii) HYPOTHETICAL: No

5

(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 69:

10

AGCCAGACTG GCTTGTCTGT AATCAGTCTA AGCTTCGGCT TGGATGCACT  
TGCAGGCTAT GCCTGCCAAC GACAATTTGG CTTGAGGGAA AAAACTAAGG  
GAAATGTGGC CCATCCGTGG GTGTTATAGT CCCTTAGAAA ATACCTTGGG  
CTGGATTGAG GTACGCAGCG AATGCTATTT GGCGAGTTGG CTGGAATAT  
15 TTTCTGAGGT GCTTT

(71) INFORMATION FOR SEQ ID NO 70:

20

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 210

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

25

(ii) MOLECULE TYPE: *Sporothrix* species specific region of 28S gene.

(iii) HYPOTHETICAL: No

30

(iv) ANTISENSE: No

## (v) SEQUENCE DESCRIPTION: SEQ ID NO: 70:

5 GACCAGACTT GCGCCYCGCG GACCACCCGG CGTTCTCGCC GGTGCACTCT  
GCGKKGCGCA GGCCAGCATC GGTTCCTCCA GGGGGACAAA GGCCGCGGGA  
ACGTAGCTCC TTCGGGAGTG TTATAGCCCG CGGCGGCATG CCCCTGGGGG  
GACCGAGGAC CGCGCTTCGG CAAGGATGCT GGCCTAATGG TCACCAGCGA  
ACCGTCTTGA

10

## (72) INFORMATION FOR SEQ ID NO 71:

## (i) SEQUENCE CHARACTERISTICS:

- 15 (A) LENGTH: 208  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

20 (ii) MOLECULE TYPE: *Scopulariopsis brevicaulis* specific region of 28S  
gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

## 25 (v) SEQUENCE DESCRIPTION: SEQ ID NO: 71:

30 GACCAGACTT GCGCCCGTCG GATCAACCGT CGCTTGCGGC GGCGCACTCC  
GGCGGGCTCA GGCCAGCATC AGTTCGTCCG GGGGGAGAAA GGCGGCGGGA  
ATGTGGCTCT TCGGAGTGTT ATAGCCCGCC GTGTAATACC CTCGGGTGGA  
CTGAGGACCG CGCGTATGCA AGGATGCTGG CGTAATGGTC GTCAGCGACC  
CGTCTTGA

## (73) INFORMATION FOR SEQ ID NO 72:

## (i) SEQUENCE CHARACTERISTICS:

- 5 (A) LENGTH: 210  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single  
(D) TOPOLOGY: linear

10 (ii) MOLECULE TYPE: *Scopulariopsis brumptii* specific region of 28S gene.

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

15

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 72:

20 GACCAGACTC GCGCCCGTCG GATCAGCCGT CGCTCGTCGG CGGCGCACTC  
CGGCGGGGCTC GGGCCAGCAT CAGTTCGCCT CGGGGGGAGA AAGGCGGCGG  
GAATGTGGCT CTACGGAGTG TTATAGCCCG CCGCGTAATA CCCCCGGGCG  
GACTGAGGAC CGCGCGTATG CAAGGATGCT GCGGTAATGG TCGTCAGCGA  
CCCGTCTTGA

25

## (74) INFORMATION FOR SEQ ID NO 73:

## (i) SEQUENCE CHARACTERISTICS:

- 30 (A) LENGTH: 214  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Saccharomyces cerevisiae* specific region of 28S gene.

5

(iii) HYPOTHETICAL: No

(iv) ANTISENSE: No

10

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 73:

15

GATCAGACAT GGTGTTTTGT GCCCTCTGCT CCTTGTGGGT AGGGGAATCT  
CGCATTTCAC TGGGCCAGCA TCAGTTTTGG TGGCAGGATA AATCCATAGG  
AATGTAGCTT GCCTCGGTAA GTATTATAGC CTGTGGGAAT ACTGCCAGCT  
GGGACTGAGG ACTGCGACGT AAGTCAAGGA TGCTGGCATA ATGGTTATAT  
GCCGCCCCGTC TTGA

20 (75) INFORMATION FOR SEQ ID NO 74:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 236

(B) TYPE: nucleic acid

25

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

(ii) MOLECULE TYPE: *Trichosporon beigelii* specific region of 28S gene.

30

(iii) HYPOTHETICAL: No



(iv) ANTISENSE: No

(v) SEQUENCE DESCRIPTION: SEQ ID NO: 74:

5

AGTCAGTCGT GTTCTTTGGA TTCAGCCAGT TCTGCTGGTC TACTTCCTTG  
GAACGGGTCA ACATCAGTTT TGTCCGGTGG ATAAAGGTAG TAGGAATGTG  
ACTTCTCCGG AAGTGTTATA GCCTATTATC ACATACACTG GGTGAGACTG  
AGGACTGCAG CTCGCCTTTA TGGCCGGCCT TCGGGCACGT TCGAGCTTAG  
10 GATGTTGACA TAATGGCTTT AAACGACCCG TCTTGA

We claim:

1. An oligonucleotide probe for the 28S subunit of fungi which is able to identify one species selected from the group consisting of *Acremonium* sp., *Aspergillus clavatus*,  
5 *Aspergillus flavus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*,  
*Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*,  
*Aspergillus ustus*, *Beauveria* sp., *Bipolaris* sp., *Blastoschizomyces* sp., *Blastomyces dermatitidis*, *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*,  
10 *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* var *gattii* serotype B, *Cryptococcus neoformans* serotype A, *Cryptococcus laurentii*,  
*Cryptococcus terreus*, *Curvularia* sp., *Fusarium* sp., *Filobasidium capsuligenum*,  
*Filobasidiella* (*Cryptococcus*) *neoformans* var *bacillispora* serotype C, *Filobasidiella* (*Cryptococcus*) *neoformans* var *neoformans* serotype D, *Filobasidium uniguttulatum*,  
15 *Geotrichum* sp., *Histoplasma capsulatum*, *Malbranchea* sp., *Mucor* sp., *Paecilomyces* sp.,  
*Penicillium species*, *Pseudallescheria boydii*, *Rhizopus* sp., *Sporothrix schenckii*,  
*Scopulariopsis brevicaulis*, *Scopulariopsis brumpti*, *Saccharomyces cerevisiae*, and  
*Trichosporon beigelii*.
- 20 2. An oligonucleotide probe of claim 1 comprising all or part of any one of the sequences (SEQ ID NO: 3) through (SEQ ID NO: 74), or any functional equivalent thereof, which is able to identify the corresponding species of fungus.
3. An oligonucleotide probe of claim 1 comprising all or part of any one of the sequences  
25 (SEQ ID NO: 3) through (SEQ ID NO: 74), or any functional equivalent thereof, which is able to identify any one or more of said species of fungus.
4. A method of determining whether one or more fungal species selected from a group of fungi is present in a sample comprising the following steps:  
30 a. extracting the nucleic acid material from the fungi contained in said sample,

- b. adding two known primers, (SEQ ID NO 1) and (SEQ ID NO 2), or the functional equivalent thereof, bracketing the areas of interest on the 28S rDNA or rRNA present in said group of interest,
  - c. amplifying the sequence between said primers, and
  - 5 d. using one or more third labeled probes to determine which of said fungi is present, wherein said third probes are selected from the group consisting of (SEQ ID NO 3) through (SEQ ID NO 74), any portion thereof and functional equivalents thereof.
- 10 5. A method of claim 4 in which said amplifying procedure is the polymerase chain reaction.
6. A method of claim 4 which, following said amplification, comprises the following step:
- d. using two or more third probes to determine which of said fungi is present, one of
- 15 said third probes being attached to a moiety which allows separation of said probe and one or more third probes connected to labeled moieties, wherein said third probes are selected from the group consisting of (SEQ ID NO 3) through (SEQ ID NO 74), any portion thereof and functional equivalents thereof.
- 20 7. A method of claim 4 which excludes said amplification step.
8. A method of claim 4 wherein said fungal species is selected from the group consisting of *Acremonium* sp., *Aspergillus clavatus*, *Aspergillus flavus*, *Aspergillus fumigatus*, *Aspergillus glaucus*, *Aspergillus nidulans*, *Aspergillus niger*, *Aspergillus ochraceus*, *Aspergillus terreus*, *Aspergillus unguis*, *Aspergillus ustus*, *Beauveria* sp., *Bipolaris* sp., *Blastoschizomyces* sp., *Blastomyces dermatitidis*, *Candida albicans*, *Candida glabrata*, *Candida guilliermondii*, *Candida kefyr*, *Candida krusei*, *Candida lusitanae*, *Candida parapsilosis*, *Candida tropicalis*, *Chrysosporium* sp., *Cladosporium* sp., *Coccidioides immitis*, *Cryptococcus neoformans* var *gattii* serotype B, *Cryptococcus*
- 25 *neoformans* serotype A, *Cryptococcus laurentii*, *Cryptococcus terreus*, *Curvularia* sp.,
- 30

*Fusarium* sp., *Filobasidium capsuligenum*, *Filobasidiella* (Cryptococcus) *neoformans* var *bacillispora* serotype C, *Filobasidiella* (Cryptococcus) *neoformans* var *neoformans* serotype D, *Filobasidium uniguttulatum*, *Geotrichum* sp., *Histoplasma capsulatum*, *Malbranchea* sp., *Mucor* sp., *Paecilomyces* sp., *Penicillium* species, *Pseudallescheria*  
5 *boydii*, *Rhizopus* sp., *Sporothrix schenkii*, *Scopulariopsis brevicaulis*, *Scopulariopsis brumpti*, *Saccharomyces cerevisiae*, and *Trichosporon beigelii*.

9. A method of claim 4 wherein more than one third probe is used, each said third probe connected to a different signal moiety or moiety which allows separation of said  
10 third probe.

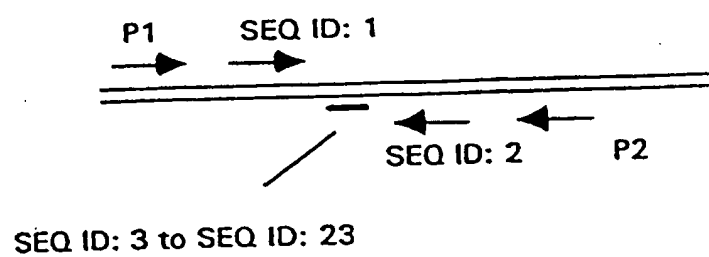


FIGURE 1

Figure 2A

1 70

```

(Rhizo2) AGCCAGACTG GCTTGTCTGT AATCAATCTA GGCTTCG.GC CTGGATGCAC TTGCAGGCTA ..TGCTGOC
(Rhizo3) AGCCAGACTG GCTTGTCTGT AATCAGTCTA AGCTTCG.GC TTGGATGCAC TTGCAGGCTA ..TGCTGOC
(Rhizo1) AGCCAGACTG GCTTGTCTGT AATCAATCTA GGTTTGTGTC CTGGATGCAC TTGCAGACTA TTTCGCTGOC
(Macor_) AGCCAGACTG GTTTGACTGT AATCAACCTA GAATTCGTTT .TGGGTGCAC TTGCAGTCTA ..TACCTGOC
(C_Terr) AGTCAGTCAT GTCTATTGGA CTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(F_Caps) AGTCAGTCAT GTCTATTGGA CTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(F_Unig) AGTCAGTCAT GTCTATTGGA CTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Neob) AGTCAGTCAT GTCTATTGGA TTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(F_Neoc) AGTCAGTCAT GTCTATTGGA TTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(F_Neod) AGTCAGTCAT GTCTATTGGA TTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Neef) AGTCAGTCAT GTCTATTGGA TTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(T_Beig) AGTCAGTCAT GTCTATTGGA TTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Laur) AGTCAGTCAT GTCTATTGGA CTCAGGCGGT TCT.....G OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Beauve) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOC. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Fusari) GACAGACTT GGGCTTGGT AATCATCTG GGTC.TCY. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Acremo) GACAGACTT GGGCTTGGT AATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Paecil) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(P_Boyd) GACAGACTT GGGCTTGGT AATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(S_Brum) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(S_Brev) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Sporot) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(B_Derm) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(H_Caps) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Nidu) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Ungu) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Ustu) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Clav) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Fum) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Flav) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Ochr) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Nige) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Tarr) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(A_Glan) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Penici) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Lumi) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Bipola) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Chitv) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Chryso) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Clados) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Malbra) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Para) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Trop) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Alb) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Gril) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Glab) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(S_Cere) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Farf) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Geotri) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Lusi) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(C_Krus) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.
(Blasch) GACAGACTT GGGCTTGGT GATCATCGG GGTC.TOG. OGGGTGACT TOCTTAGAT GGGGTCAAC.

```

71 140

(Shizo2)	AACGACAATT	TGACTTGAGG	GAJAJAACTA	GGGGAAATGT	GGCC.....	CACTTGTTGGG	TGTTATAGTC
(Rhizo3)	AACGACAATT	TGGCTTGAGG	GAJAJAACTA	AGGGAAATGT	GGCC.....	CATCOGTGGG	TGTTATAGTC
(Rhizo1)	AACGACAATT	TTTTTTGAGT	GTAAJAACTA	TTGGAAATGT	GGCCAAATAT	TATTIATTTG	TGTTATAGTC
(Mucor)	AACACAGTT	TGATTTGAGG	GAJAJAACTA	GTAGGAATGT	AGCC.....	.....TCTOGA	GGTGTATAG
(C_Terr)	.ATCAGTTTT	.GATCGCTGG	AAJAGGGCAG	GAGGAATGTA	GCATC.TOG	GGTGAACCTA	TAGCCTTCTG
(F_Caps)	.ATCAGTTTT	.GACCGTTGG	ATAAAGGCAG	GAAGAATGTA	GCATC.TOG	GGTGAACCTA	TAGCCTTCTG
(F_Unig)	.ATCAGTTTT	.GATCGCTGG	ATAAAGGCAG	GAGGAATGTA	GCACCC.OOG	GGTGAACCTA	TAGCCTTCTG
(C_Neob)	.ATCAGTTCT	.GATCGGTGG	ATAAGGGCTG	GAGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(F_Neoc)	.ATCAGTTCT	.GATCGGTGG	ATAAGGGCTG	GAGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(F_Neod)	.ATCAGTTCT	.GATCGGTGG	ATAAGGGCTG	GAGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Neof)	.ATCAGTTCT	.GATCGGTGG	ATAAGGGCTG	GGGAATGTA	GCACCTCTOG	GAGTGTGTTA	TAGCCTTCTG
(T_Beig)	.ATCAGTTTT	.GTCCGGTGG	ATAAGGGCTG	TAGGAATGTA	.....ACTCTOC	GGAGTGTGTA	TAGCCTTCTG
(C_Laur)	.ATCAGTTTT	.GTCCGAAGG	ATAAGGGCTG	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Beanve)	.ATCAGTTCT	.OCC.T.GGGGG	ACAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Fusari)	.ATCAGTTTT	.CSOC.GGGGG	ACAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Acremo)	.ATCAGTTCT	.OCC.GGGGG	ACAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Paecil)	.ATCAGTTCT	.OCC.GGGGG	ACAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(P_Boyd)	.ATCAGTTCT	.CTGACGGGG	AGAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(S_Brum)	.ATCAGTTCT	.OCTCGGGGG	AGAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(S_Brev)	.ATCAGTTCT	.TOCGGGGG	AGAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Sporot)	.ATCGTTCT	.C.CCAGGGGG	ACAAAGGCCT	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(B_Dema)	.GTGGTTTC	.GACGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(H_Caps)	.GTGGTTTC	.GACGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Nidu)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Ungu)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Usta)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Clav)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Fum)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Elav)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Ochr)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Nige)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Ferr)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(A_Glan)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Penici)	.GTGGTTTC	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Imul)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Bipola)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Curvul)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Chryso)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Clados)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Malbra)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Baza)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Trop)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Alib)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Gril)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Glab)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(S_Cere)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Kefy)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Geotri)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Lind)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(C_Krus)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG
(Blasch)	.ATCAGTTCT	.GGGGGGGG	TCAAAGGCTC	GGGAATGTA	GCACCTCTOG	GGGTGTGTTA	TAGCCTTCTG

141

210

(Rhizo2) CCTTAGAANA TACCTTGGGT TGGATTGAGG AACGCAGCGA ATG..... CTTATTG  
 (Rhizo3) CCTTAGAANA TACCTTGGGC TGGATTGAGG TAGGCAGCGA ATG..... CTATTTG  
 (Rhizo1) CCTTAGAANA TACCTTGAAT TGGATTGAGG AACGCAGCGA ATGCTTCTCT TTnGAGGCAA AGTCTTTTAT  
 (Hucor\_) OCTACTATCA TACTCTGGAT TGGACTGAGG AAGGCAGCGA ATGCGWTAG GCRAGATTGC TGGGTGCTTT  
 (C\_Terr) TGGTATACAG TGGTITGGGAC TGAGGAAGGC AGCATGOCIT TATGGGCGGG GTTGGGCGAC GTACATGCTT  
 (F\_Caps) TCACATACAA TGGTITGGGAC TGAGGAAGGC AGCATGOCIT TATGGGCGGG ATTGGTGCAC GTACATGCTT  
 (F\_Unig) TCACATACAG TGGTITGGGAC TGAGGAAGGC AGCATGOCIT TATGGGCGGG ATTGGTGCAC GTACATGCTT  
 (C\_Neob) TGGCATACAC TGGTITGGGAC TGAGGAAGGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (F\_Neoc) TGGCATACAC TGGTITGGGAC TGAGGAAGGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (F\_Neod) TGGCATACAC TGGTITGGGAC TGAGGAAGGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (C\_Neof) TGGCATACAC TGGTITGGGAC TGAGGAAGGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (T\_Beig) TCACATACAC TGGGTGAGAC TGAGGACTGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (C\_Laur) TGGCATAGC GGGATGAGAC TGAGGACTGC AGCTGOCIT TATGGGCGGG GTTGGGCGAC GTTGGAGCIT  
 (Beauve) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CATTOGCA  
 (Fusari) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CHTCTGCA  
 (Acremo) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... C.TCTGCA  
 (Paecil) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... C.TGCGCA  
 (P\_Boyd) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CATCTGCA  
 (S\_Brum) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CGTATGCA  
 (S\_Brev) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CGTATGCA  
 (Sporot) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (B\_Derm) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (H\_Caps) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (A\_Nido) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (A\_Ungu) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (A\_Ustu) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCA  
 (A\_Clav) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCT  
 (A\_Fum) .....G TGGTATAGOC OGTTGGGTAA TACG.CTGG GGGGACTGAG GTTGGG... CITTOGCT  
 (A\_Flav) .....A CCTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (A\_Ochr) .....G CCTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (A\_Nige) .....A CCTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (A\_Terr) .....G CCTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (A\_Glan) .....G CCTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (Penici) .....G TGGTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCT  
 (C\_Imml) .....G TGGTATAGOC GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (Bipola) .....G CCTATAG.G GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CATCTGCT  
 (Curvul) .....G CCTATAG.G GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CATCTGCT  
 (Chryso) .....T GCTATAG.C TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCT  
 (Clados) .....T TATA.G CCTATAG.G TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCT  
 (Halbra) .....G TGGTATAGOC TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (C\_Para) .....G TGGTATAGOC T.TGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CTTOG.GCT  
 (C\_Trop) .....G TGGTATAGOC T.TGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TTTAT.AOCT  
 (C\_Albi) .....G TGGTATAGOC T.TGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TTTTAAOCT  
 (C\_Gril) .....G TGGTATAGOC T.GGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TTTT.AICA  
 (C\_Glab) .....G TGGTATAGOC C.TGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CTGGTATCT  
 (S\_Cere) .....G TGGTATAGOC T.GGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TGA...TCA  
 (C\_Kefy) .....G TGGTATAGOC C.GGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TTTT...TCA  
 (Geotri) .....G TGGTATAGOC T.TGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... CITTOGCA  
 (C\_Insl) .....G TGGTATAGOC C.GGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... TCT  
 (C\_Krus) .....G TGGTATAGOC A.GGCT.GA TAGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... GTGAGGICA  
 (Blasch) .....TGAA TGGTGAAG GGGGGTGCAA TGCGGCGAGC CTGGAGCGAG GAAAGG... GTGAGGICA



	211	250
{Rhizo2}	GCGAGTTTTC CAGGAAGCT. ....TTTCT	GAGGTACTAC
{Rhizo3}	GCGAGTTGGC TGGGAATAT. ....TTTCT	GAGGTGCTTT
{Rhizo1}	TGGGATTTAC GGATCAGAC. ....TGTTG	CATTGTCACA
{Mucor_}	CGCTAATAAA TGTTAGAATT TCTGCTTGG	GTGGTGCTAA
{C_Terr}	AGG..ATGTT GACATAATGG CTTTAAAOGA	COOGTCITGA
{F_Caps}	AGG..ATGTT GACATAATGG CTTTAAAOGA	COOGTCITGA
{F_Unig}	AGG..ATGTT GACATAATGG CTTTAAAOGA	COOGTCITGA
{C_Neob}	AGG..ATGTT GACAAAATGG CTTTAAAOGA	COOGTCITGA
{F_Neoc}	AGG..ATGTT GACAAAATGG CTTTAAAOGA	COOGTCITGA
{F_Neod}	AGG..ATGTT GACAAAATGG CTTTAAAOGA	COOGTCITGA
{C_Neof}	AGG..ATGTT GACAAAATGG CTTTAAAOGA	COOGTCITGA
{T_Belg}	AGG..ATGTT GACATAATGG CTTTAAAOGA	COOGTCITGA
{C_Laur}	AGG..ATGTT GAOGTAATGG CTTTAAAOGA	COOGTCITGA
{Beauve}	AGG..ATGCT GGOGTAATGG TCATCAGTGA	COOGTCITGA
{Fusari}	AGG..ATGCT GGOGTAATGG TCATCAGTGA	COOGTCITGA
{Acremo}	AGG..ATGCT GGOGTAATGG TCATCAGTGA	COOGTCITGA
{Paecil}	AGG..ATGCT GGOGTAATGG TCATCAGTGA	COOGTCITGA
{P_Boyd}	AGG..ATGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{S_Brum}	AGG..ATGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{S_Brev}	AGG..ATGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{Sporot}	AGG..ATGCT GGOGTAATGG TCACCAAGGA	COOGTCITGA
{B_Dema}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{H_Caps}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{A_Nidu}	CGG..AOGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{A_Ungu}	CGG..AOGCT GGCTTAATGG TTGCARAAGA	COOGTCITGA
{A_Ustu}	CGG..AOGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{A_Clav}	CGG..AOGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{A_Fumi}	CGG..AOGCT GGOGTAATGG TOGTCAAGGA	COOGTCITGA
{A_Flav}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{A_Ochr}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{A_Nige}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{A_Terr}	CGG..AOGCT GGCTTAATGG TTGTCAAGGA	COOGTCITGA
{A_Glau}	CGG..AOGCT GGCTTAATGG TOGTCAAGGA	COOGTCITGA
{Penici}	CGG..AOGCT GGCTTAATGG TOGTAAA...	.....
{C_Imai}	CGG..ATGCT GGCTTAATGG TTGTCAAGGG	COOGTCITGA
{Bipola}	AGG..ATGCT GGOGTAATGG CIGTAAGGG	COOGTCITGA
{Curvul}	AGG..ATGCT GGOGTAATGG CIGTAAGGG	COOGTCITGA
{Chryso}	AGG..ATGCT GGOGTAATGG TTGTCAAGGG	COOGTCITGA
{Clados}	AGG..ATGCT GGOGTAATGG TOGTCAATGG	COOGTCITGA
{Malbra}	CGG..ATGCT GGOGTAATGG CIGTAAGGG	COOGTCITGA
{C_Para}	AGG..ATGCT GGCTTAATGA TCCTAAGTGG	COOGTCITGA
{C_Trop}	AGG..ATGCT GGCTTAATGA TCCTAAGTGG	COOGTCITGA
{C_Albi}	AGG..ATGCT GGCTTAATGA TCCTAAGTGG	COOGTCITGA
{C_Gd1}	AGG..ATGCT GGCTTAATGA TOGTCAAGGG	COOGTCITGA
{C_Glab}	AGG..ATGCT GGCTTAATGG TTGTCAAGGG	COOGTCITGA
{S_Cere}	AGG..ATGCT GGCTTAATGG TTGTCAAGGG	COOGTCITGA
{C_Kofy}	AGG..ATGCT GGOGTAATGG TTAAATGGGG	COOGTCITGA
{Geotri}	AGG..AOGCT GGCTTAATGA TTCTATAOOG	COOGTCITGA
{C_Imai}	AGG..AOGCT GGOGTAATGG TTGTCAAGGG	COOGTCITGA
{C_Kms}	OGG..ATGCT GGCTTAATGG OGCTCAAGGG	COOGTCITGA
{Blasch}	GGGCGAGGT CTCCTGAGGC TCCTTTGGCA	COOGTCITGA

# INTERNATIONAL SEARCH REPORT

Inter national Application No  
PCT/IB 96/00026

A. CLASSIFICATION OF SUBJECT MATTER  
IPC 6 C12Q1/68

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)  
IPC 6 C12Q C07K

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	JOURNAL OF MEDICAL AND VETERINARY MYCOLOGY, 32 (5). 1994. 331-341., XP002002929 LECLERC M C ET AL: "Phylogeny of dermatophytes and dimorphic fungi based on large subunit ribosomal RNA sequence comparisons" cited in the application	1
A	see the whole document	2-9
X	CURRENT GENETICS, 27 (1). 1994. 38-45., XP002002930 NEUVEGLISE C ET AL: "Identification of group-I introns in the 28s rDNA of the entomopathogenic fungus Beauveria brongniartii"	1
A	see the whole document	2-9
	-/--	

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

### \* Special categories of cited documents:

- \*A\* document defining the general state of the art which is not considered to be of particular relevance
- \*E\* earlier document but published on or after the international filing date
- \*L\* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \*O\* document referring to an oral disclosure, use, exhibition or other means
- \*P\* document published prior to the international filing date but later than the priority date claimed

- \*T\* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- \*X\* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- \*Y\* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- \*A\* document member of the same patent family

Date of the actual completion of the international search

21 May 1996

Date of mailing of the international search report

04.06.1996

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentaan 2  
NL - 2280 HV Rijswijk  
Tel. (+ 31-70) 340-2040, Tx. 31 651 epo nl,  
Fax (+ 31-70) 340-3016

Authorized officer

Gurdjian, D

## INTERNATIONAL SEARCH REPORT

Inter. Application No.  
PC1/IB 96/00026

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT		
Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	PHYTOPATHOLOGY, 84 (3). 1994. 256-259., XP002002931 MOUKHAMEDOV R ET AL: "Use of polymerase chain reaction-amplified ribosomal intergenic sequences for the diagnosis of Verticillium tricorpus"	1
A	see the whole document ---	2-9
X	CURR GENET, 12 (3). 1987. 209-214., XP002002932 CARR L G ET AL: "ORGANIZATION OF THE 5.8S 16-18S AND 23-28S RIBOSOMAL RNA GENES OF CEPHALOSPORIUM-ACREMONIUM"	1
A	see the whole document ---	2-9
X	JOURNAL OF BACTERIOLOGY, vol. 172, no. 8, 1990, pages 4238-4246, XP002002933 VILGALYS R. ET AL.: "Rapid genetic identification and mapping of enzymatically amplified ribosomal DNA from several Cryptococcus species"	1
A	see the whole document ---	2-9
A	US,A,5 352 579 (MILLIMAN CURT L) 4 October 1994 see the whole document ---	2-9
A	EP,A,0 422 872 (GENE TRAK SYSTEMS) 17 April 1991 see the whole document ---	2-9
P,X	JOURNAL OF CLINICAL MICROBIOLOGY, 33 (11). 1995. 2913-2919., XP002002934 SANDHU G S ET AL: "Molecular probes for diagnosis of fungal infections" see the whole document -----	1-3

# INTERNATIONAL SEARCH REPORT

International Application No  
PCT/IB 96/00026

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
US-A-5352579	04-10-94	NONE	
EP-A-0422872	17-04-91	AU-B- 6390490	18-04-91
		CA-A- 2025181	13-04-91
		JP-A- 3168085	19-07-91
		US-A- 5324632	28-06-94